

SCHOOL ON SYNCHROTRON RADIATION AND APPLICATIONS
In memory of J.C. Fuggle & L. Fonda

19 April - 21 May 2004

Miramare - Trieste, Italy

1561/11

Optical components for hard x-ray beamlines

Edoardo Busetto


abdus salam
international centre for theoretical physics

School on Synchrotron Radiation

27 -28 April 2004

Optical components for hard x-ray beamlines

Edoardo Busetto

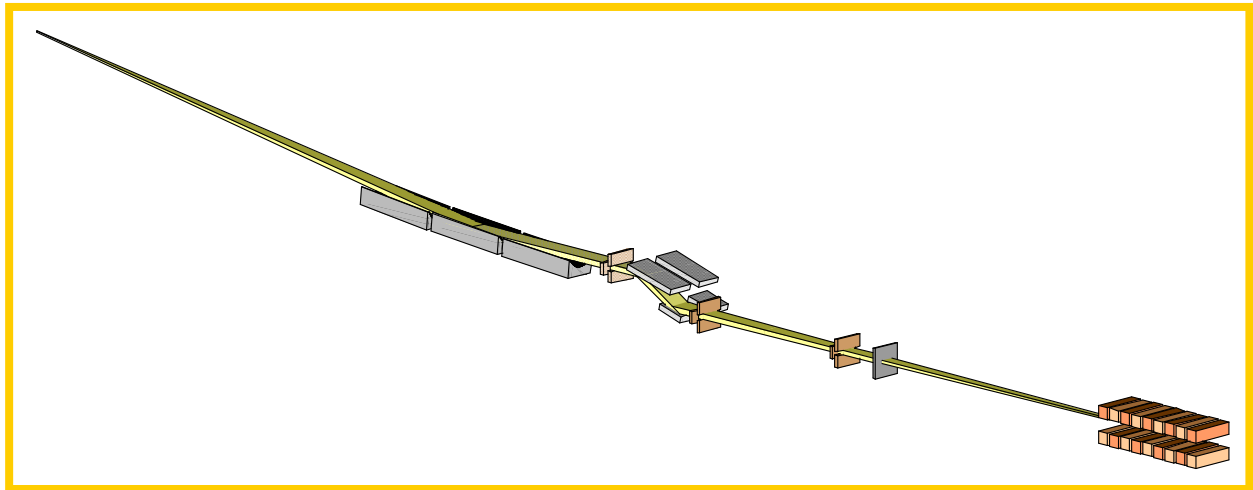

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Optical components for hard x-ray beamline


High brilliance and small electron beam emittance mean X-ray beams of high quality

“ ...The finite quality and the fundamental limits of the optical components increase the emittance of the beam.....

The main aim of the optical design consists on minimizing the inevitable beam degradation ”

Jean Susin "Design parameters for hard x-ray mirrors: the ESRF case"
OPTICAL ENGINEERING/February 1995/Vol. 34 2/361

detectors

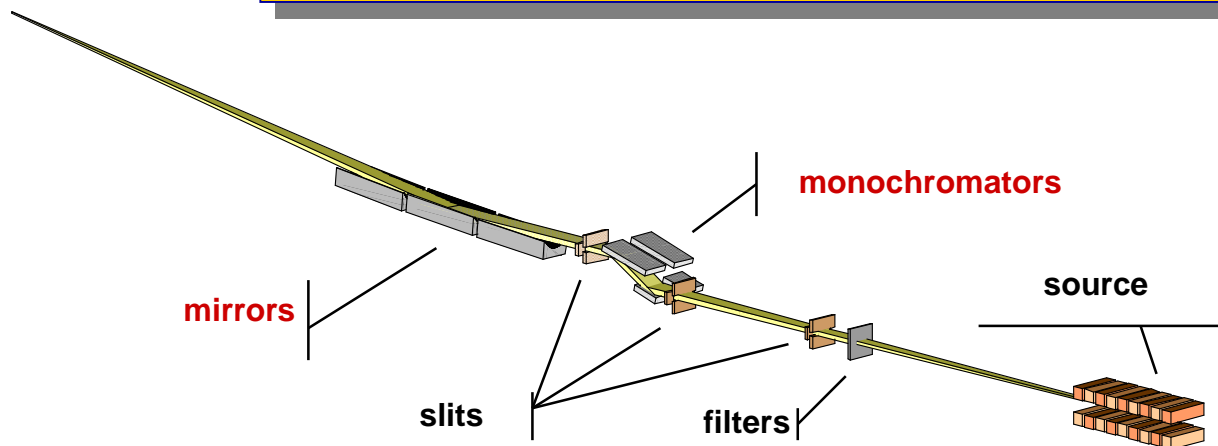

Most important optical elements:

- sources
- filters
- slits and pinholes
- mirrors
- monochromators
- detectors




Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Mirrors 1: total reflection

For x-rays the refractive index is $n = 1 - \delta$
with $0 < \delta \ll 1$, therefore is $0 < n < 1$

If we consider ϑ as the angle that the incoming radiation does with the mirror surface (*grazing angle*), the photons will be totally reflected if $\vartheta < \vartheta_c$

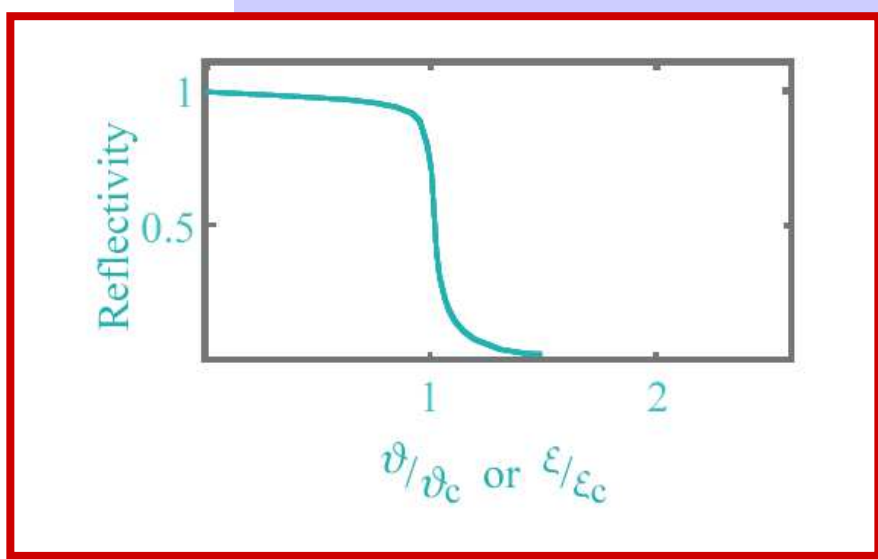


$$\vartheta < \vartheta_c$$




$$\vartheta = \vartheta_c$$



$$\vartheta > \vartheta_c$$

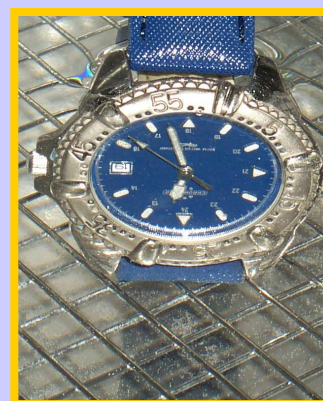

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

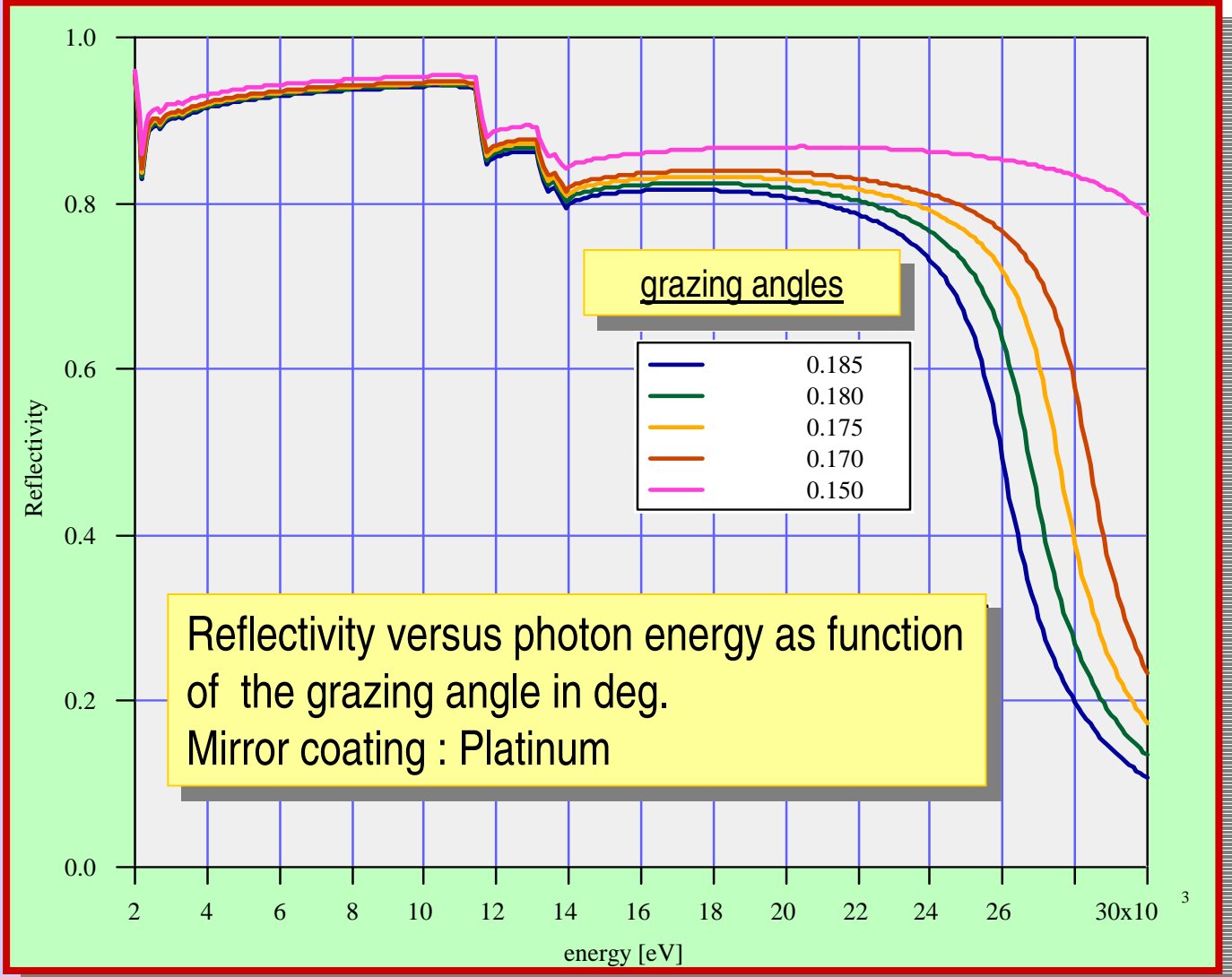

Reflectivity versus photon energy as function of the grazing angle in deg.
 Mirror coating : Platinum




Sincrotrone Trieste S.c.p.A.
 Hard-X Ray Optics
 Laboratory

Mirrors 2: focussing

In the ideal mirror device all rays from one particular **point** are reflected and focused into another **point** according to $1/q + 1/p = 1/f$


ELLIPSOIDAL SURFACE

best approximation circle

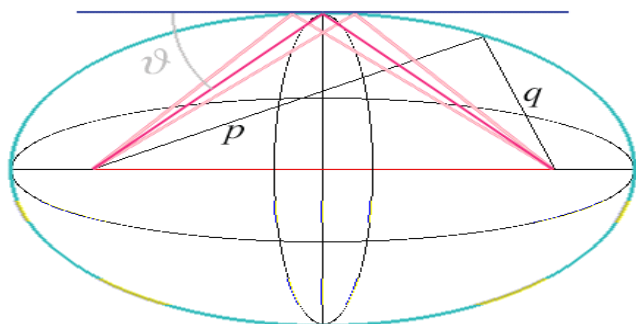

$$\rho_{\text{tangential}} = 2f / \sin \vartheta$$

Tangential focussing

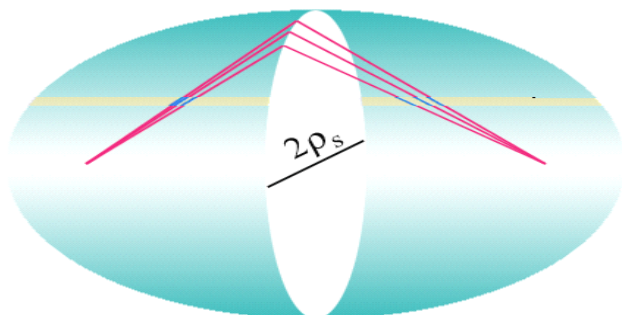




$$\rho_{\text{sagittal}} = 2f \sin \vartheta$$


Sagittal focussing

- * the bending magnet case
- * the extended source case



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

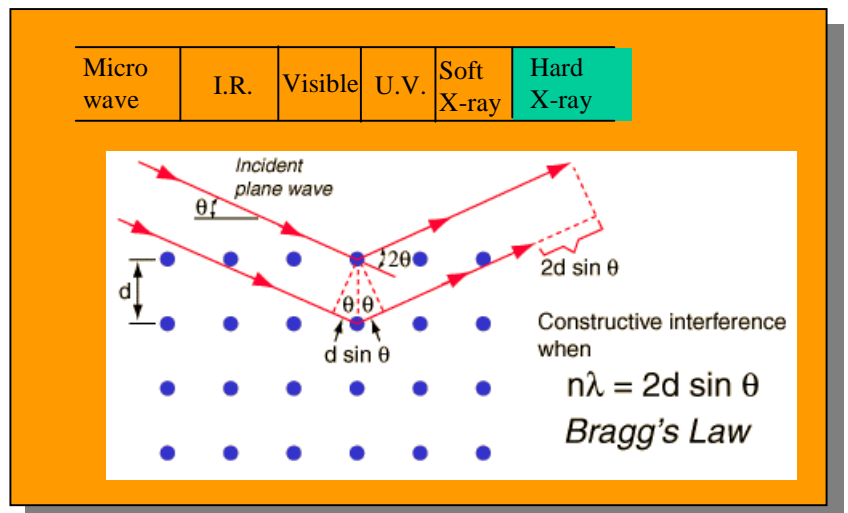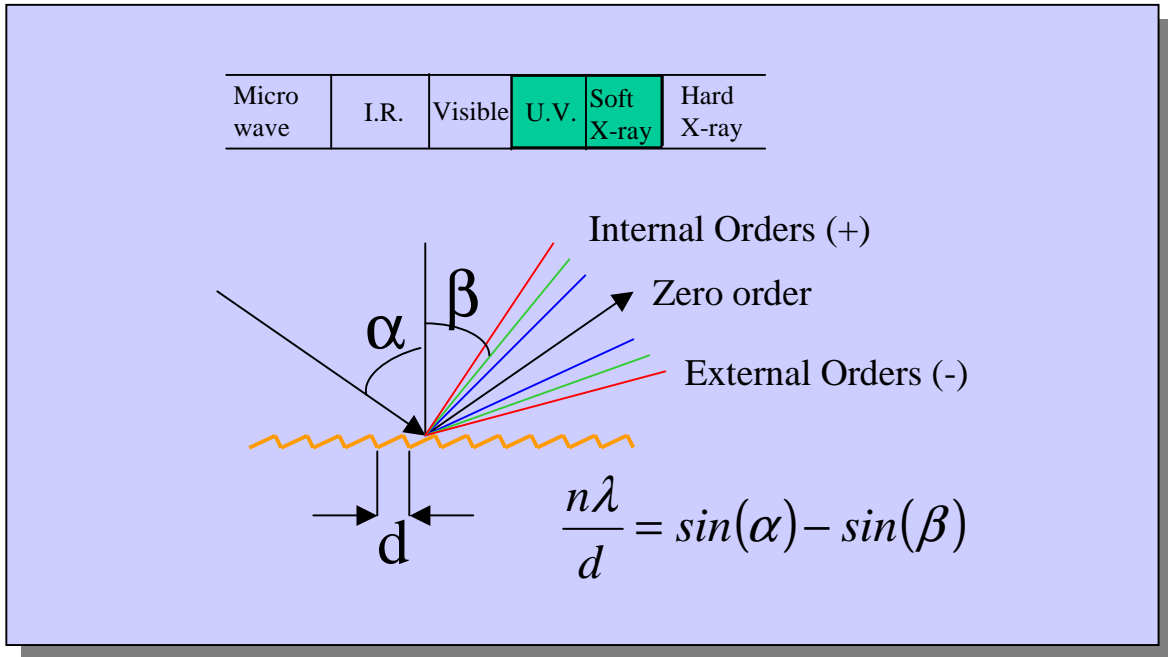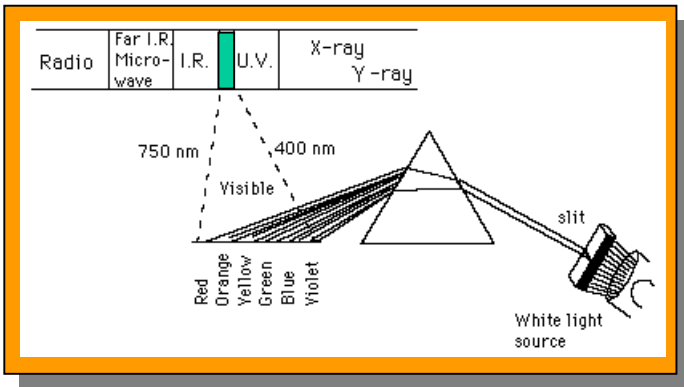

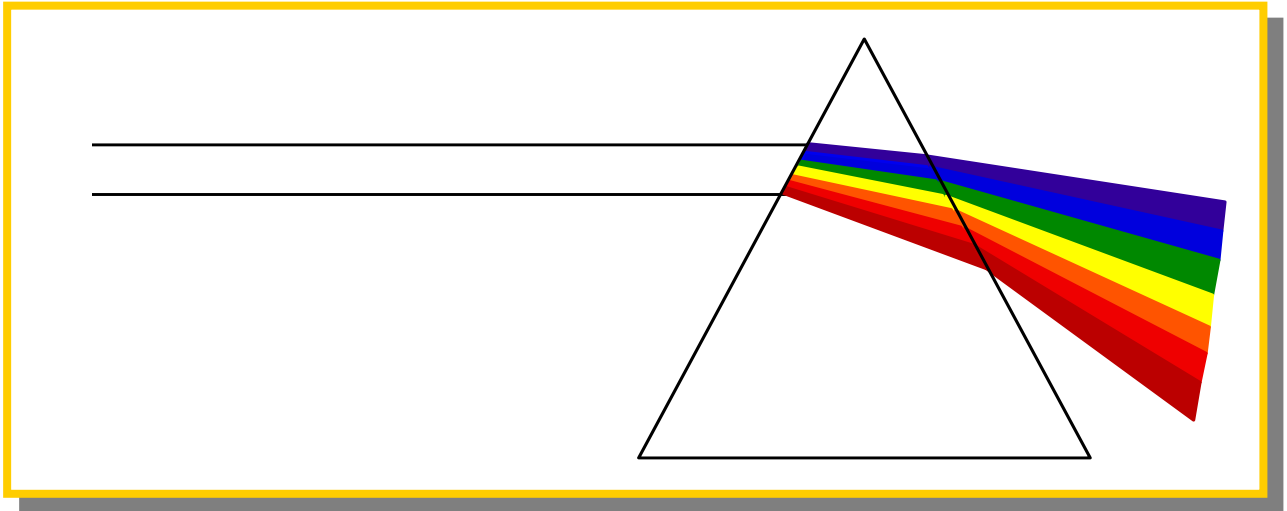


the optical prism is used to separate the components of the white visible light

sampling the out coming light with a slit it is possible to select a part of the spectrum with a spectral purity which depends on the distance and the slits aperture.



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

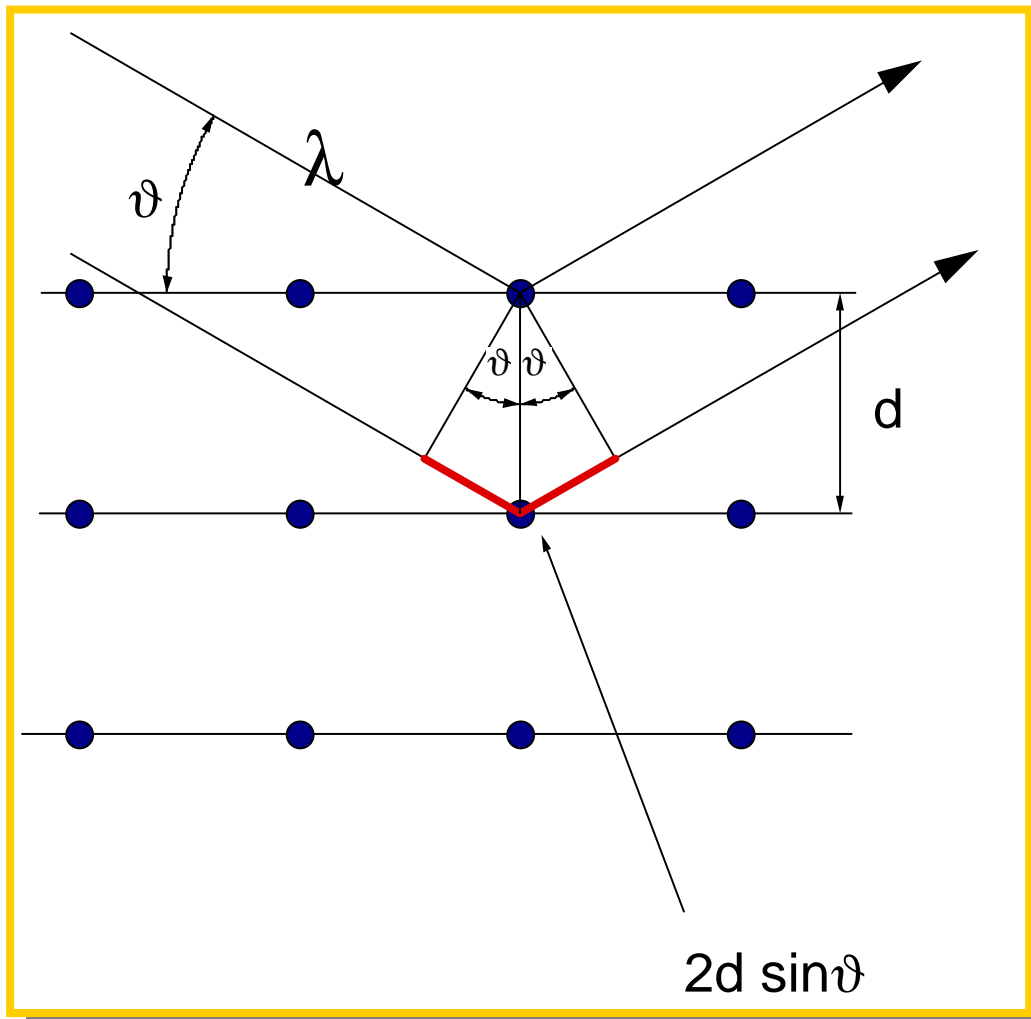

the Bragg's law

Radiation of wavelength λ is reflected by the lattice plane.

The outgoing waves interfere. The interference is constructive only if the difference of optical path is a multiple of λ :

$$2d \sin \vartheta = n \lambda$$



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

$$2d\sin\vartheta = n\lambda$$

from the Bragg law

$$\sin\vartheta = 1 \Rightarrow \lambda_{\max}$$

therefore

$$\lambda_{\max} = 2d$$

and the Bragg angle is 90°



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

important properties for the x-ray monochromators

- *ENERGY RESOLUTION*

$$\frac{\Delta\lambda}{\lambda} = \frac{\Delta E}{E} = \Delta\vartheta \cot g(\vartheta_B)$$

$\Delta\vartheta$ has two contribution :


$\Delta\vartheta_{\text{beam}}$ - beam angular spread (optics)

ω_{crystal} - intrinsic reflection width of the
monochromator




Case of $\Delta\vartheta_{beam} \gg \omega_{crystal}$
white beam with divergence in the plane of scattering

The crystal accepts all the rays with
 $\vartheta_{min} \leq \vartheta_B \leq \vartheta_{max}$

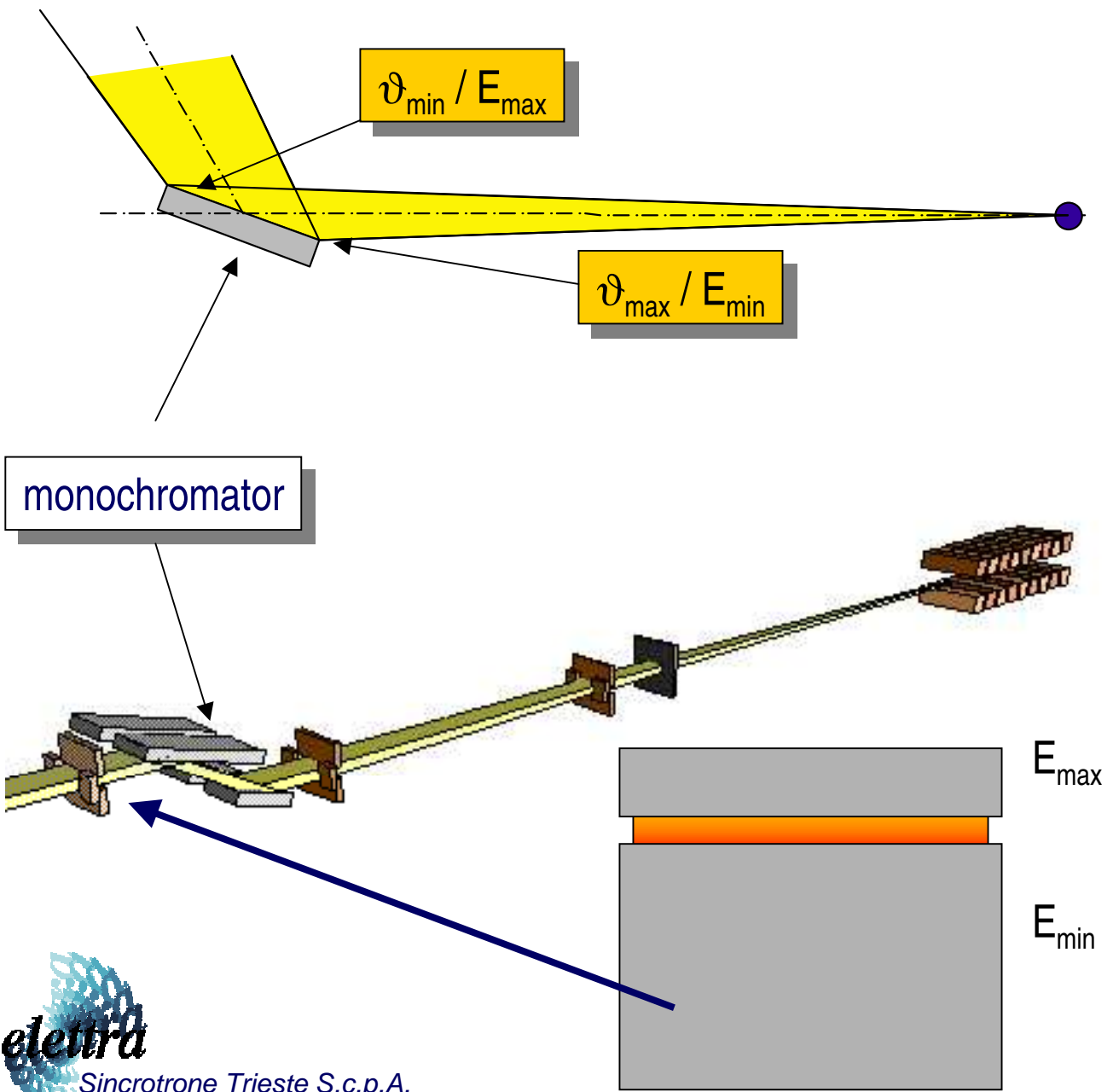

Case of $\Delta\vartheta_{beam} \ll \omega_{crystal}$
white beam parallel in the plane of scattering

parabola

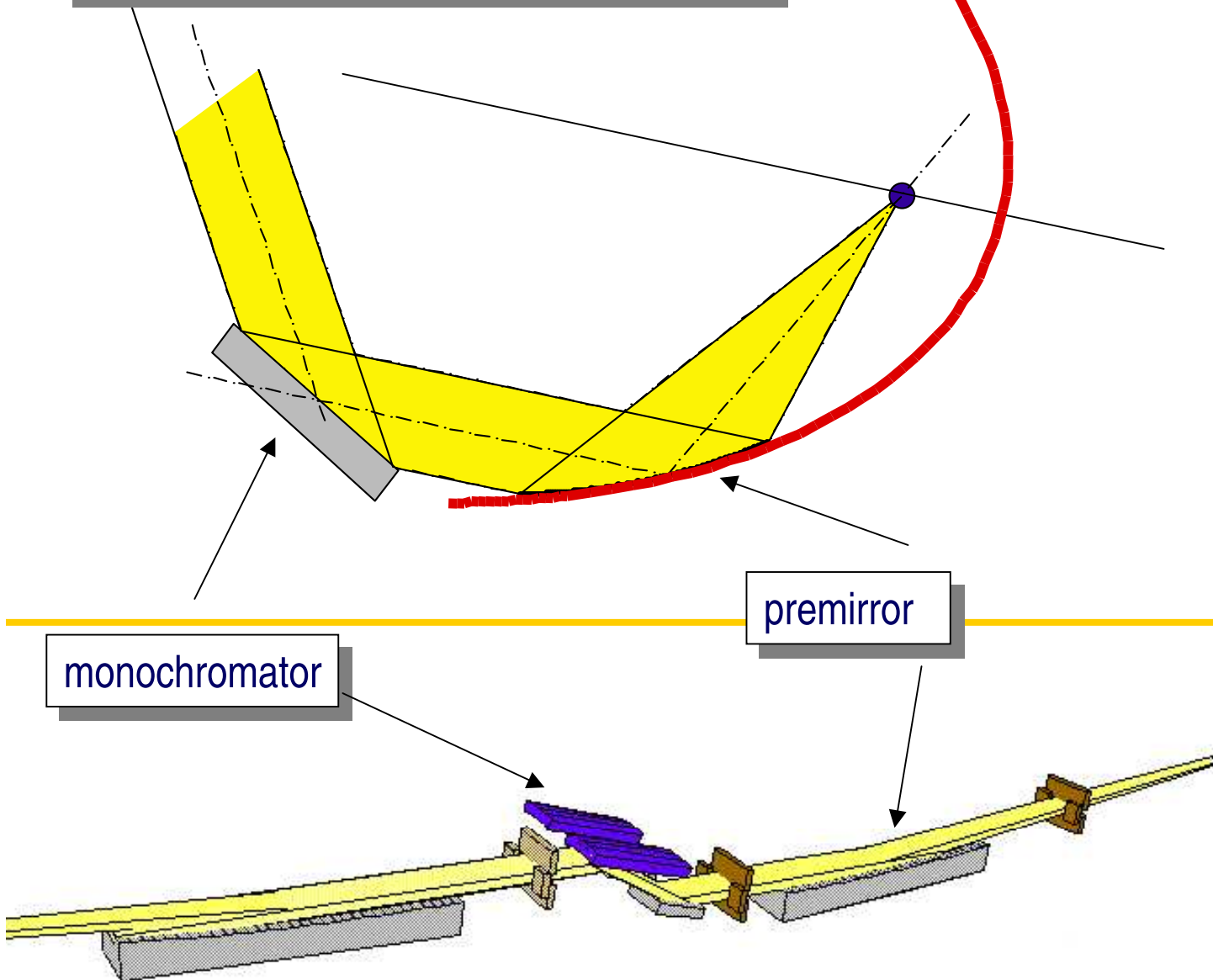


The energy bandwidth is determined by the derivative of the Bragg's law

$$\Delta E = \omega_{crystal} \cotg(\vartheta_B) E$$





Sincrotrone Trieste S.c.p.A.
 Hard-X Ray Optics
 Laboratory


Dumond diagrams

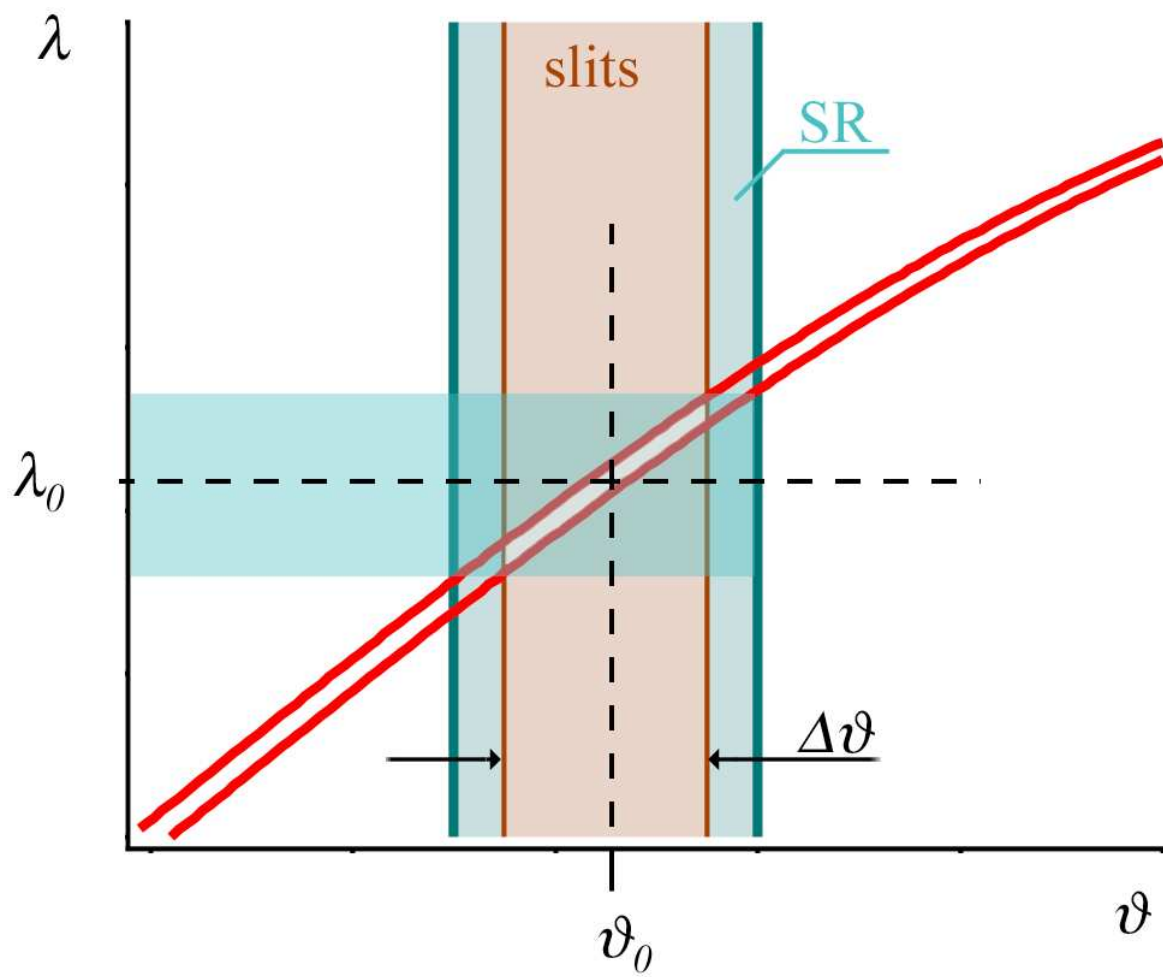

Second crystal in
non dispersive configuration

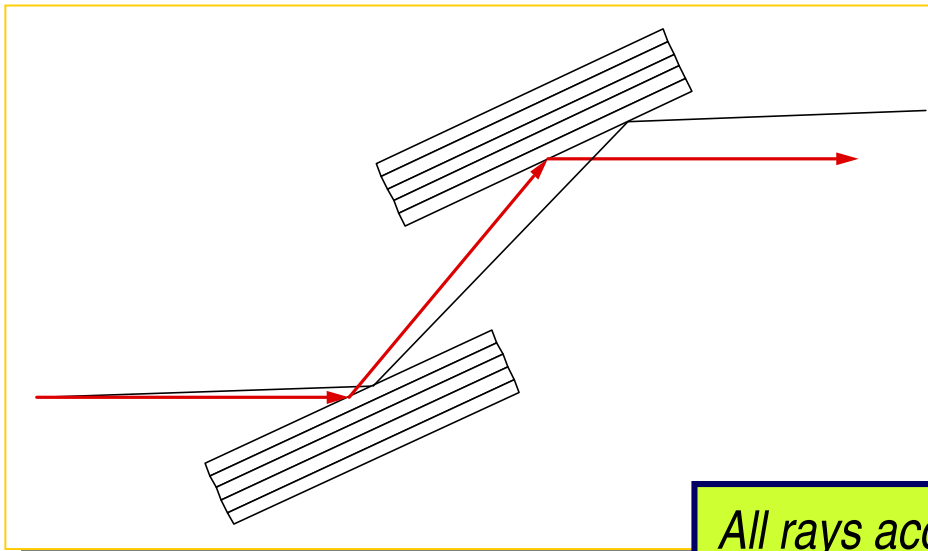
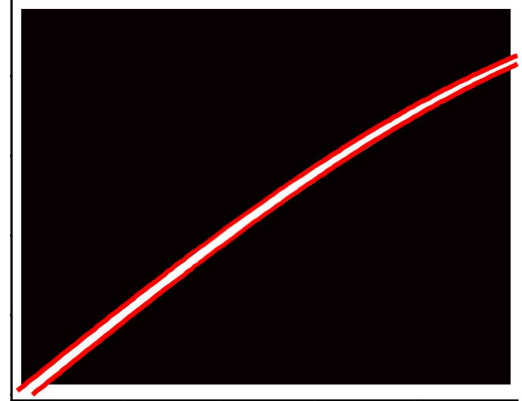


*All rays accepted
by the first crystal
are accepted also
at the second.*

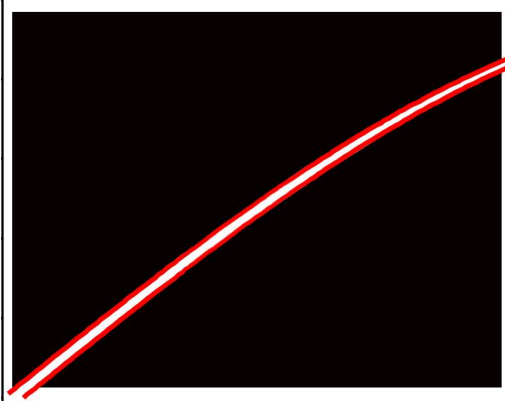


Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Second crystal in **dispersive** configuration

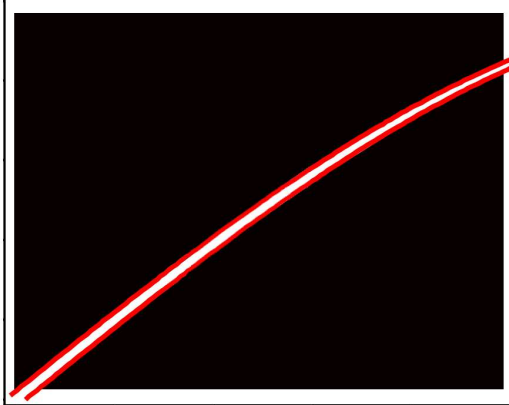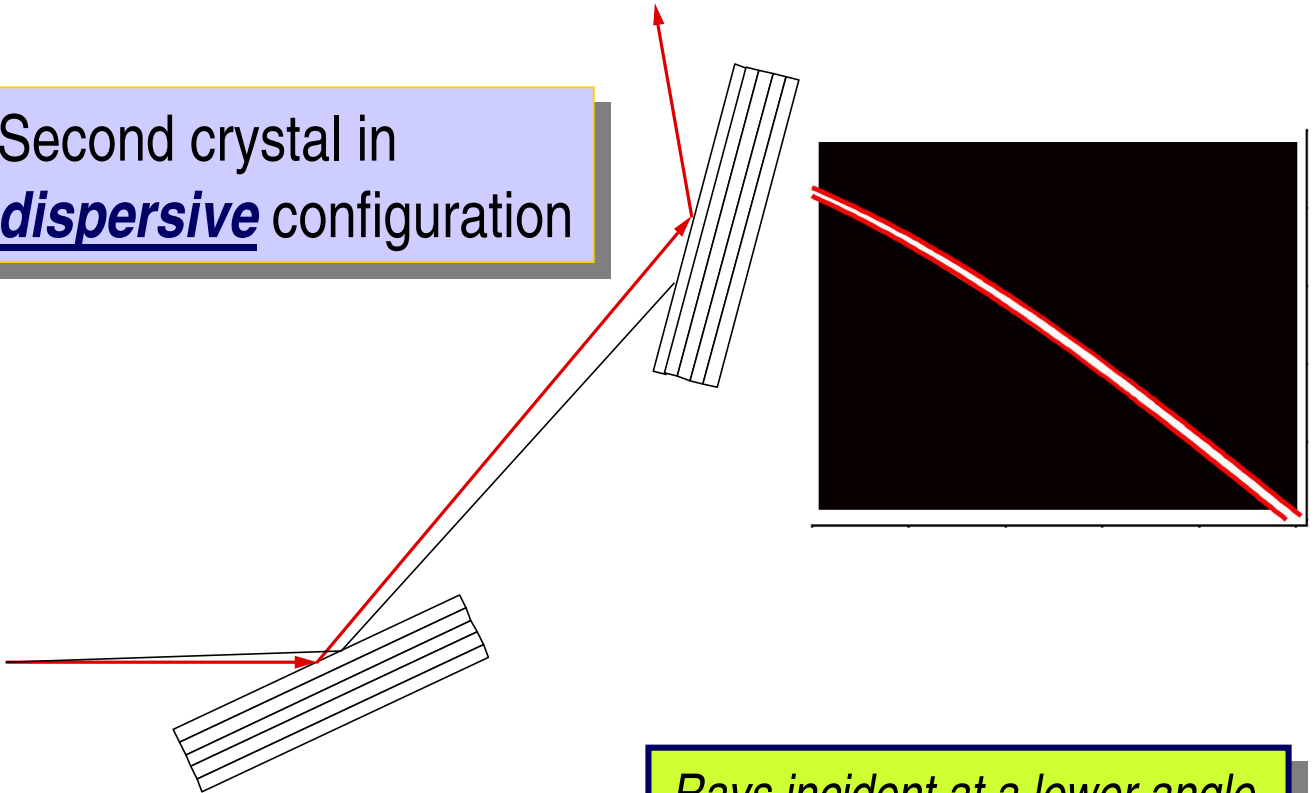

Rays incident at a lower angle than the central ray on the first crystal are incident at a higher angle on the second crystal.

energy resolution ↑
intensity of the reflection ↓



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

two models for the x-ray diffraction in single crystals

kinematical model

apply this model for:

- **thin perfect crystals**
- **distorted or mosaic crystals**

according with Darwin model (1922) the mosaic crystal is defined by two general conditions:

- crystallites have to be **misoriented** more than the Darwin width of the perfect crystal (loss of the phase condition)
- their dimensions have to be smaller than the **extinction length** of the considered radiation (no second interaction)

dynamical model

apply this model for:

- **thick and perfect crystal**

a) we can't longer consider single interaction. (extinction length)

b) we can't neglect, as well as in the kinematical model, the effect of the radiation absorption



INTENSITY OF THE REFLECTION

- reflectivity or peak reflectivity
- integral reflecting power

The Darwin curve


$$\omega_s = \frac{2}{\sin 2\vartheta_B} \frac{r_e \lambda^2}{\pi V} C |F_{hr}| e^{-M}$$

n order of the reflection
 λ_1 wavelength of the fundamental
 $e^{-M(n)}$ temperature factor
 V volume of the unit cell
 ϑ_B Bragg angle
 R_e radius of the electron e^2/mc^2

F_{hr} real part of the structure factor related to the diffracted direction $\mathbf{h}(h,k,l)$

predicted by the dynamical theory

absorption effects

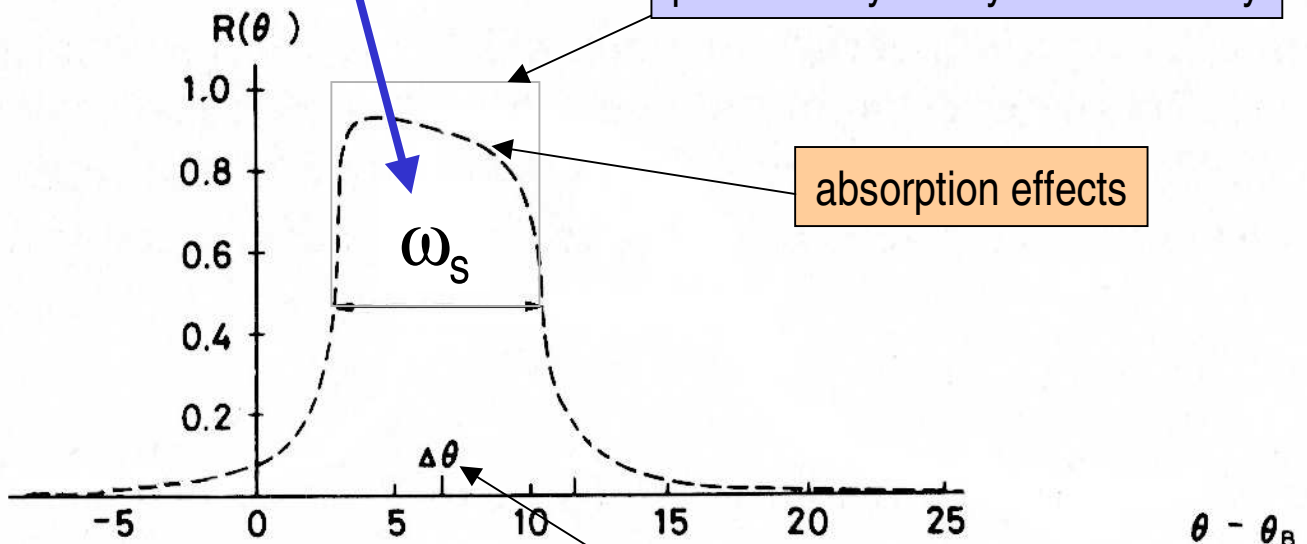

angular shift due to the refractive effect



the **b** parameter defined as :

$$b = \frac{\sin(\alpha - \vartheta_B)}{\sin(\alpha + \vartheta_B)}$$

α is the angle between the Bragg plane and the crystal surface

T. Matsushita and H. Hashizume X-Ray Monochromators Handbook on Synchrotron Radiation, Vol. 1

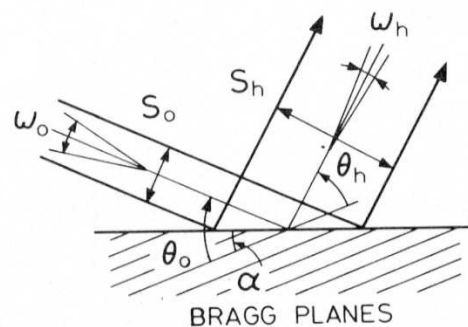


Fig. 3. Geometry of X-ray reflection by a perfect single crystal. θ_0 : incidence angle; θ_h : reflection angle. For a non-zero asymmetry angle α ($0 < |\alpha| < \theta_B$), the angular width ω_0 for acceptance is not equal to the angular width ω_h for emergence. The figure is drawn for $b < 1.0$, where $\omega_0 > \omega_s > \omega_h$. Note also the change of beam cross sections, S_0 and S_h .

$$\omega_0 = \frac{\omega_s}{\sqrt{b}}$$

the angular acceptance as function of the intrinsic width and the **b** parameter:



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

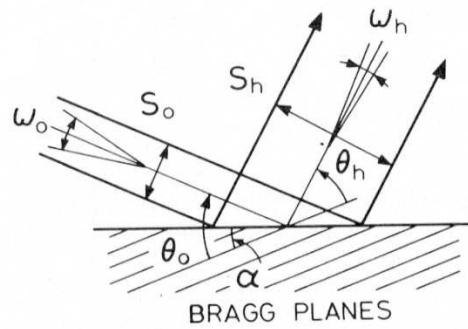


Fig. 3. Geometry of X-ray reflection by a perfect single crystal. θ_0 : incidence angle; θ_h : reflection angle. For a non-zero asymmetry angle α ($0 < |\alpha| < \theta_B$), the angular width ω_0 for acceptance is not equal to the angular width ω_h for emergence. The figure is drawn for $b < 1.0$, where $\omega_0 > \omega_s > \omega_h$. Note also the change of beam cross sections, S_0 and S_h .

Bragg reflection width in case of asymmetric cut crystal is defined by:

$$\omega_h = \omega_s \sqrt{b}$$

$$\omega_h = b \omega_0$$

the angular acceptance as function of the Bragg reflection width



also for the beams sections

$$S_h = \frac{S_0}{b}$$

combining the two formulas we have the well known Liouville's theorem

$$\omega_h S_h = \omega_0 S_0$$



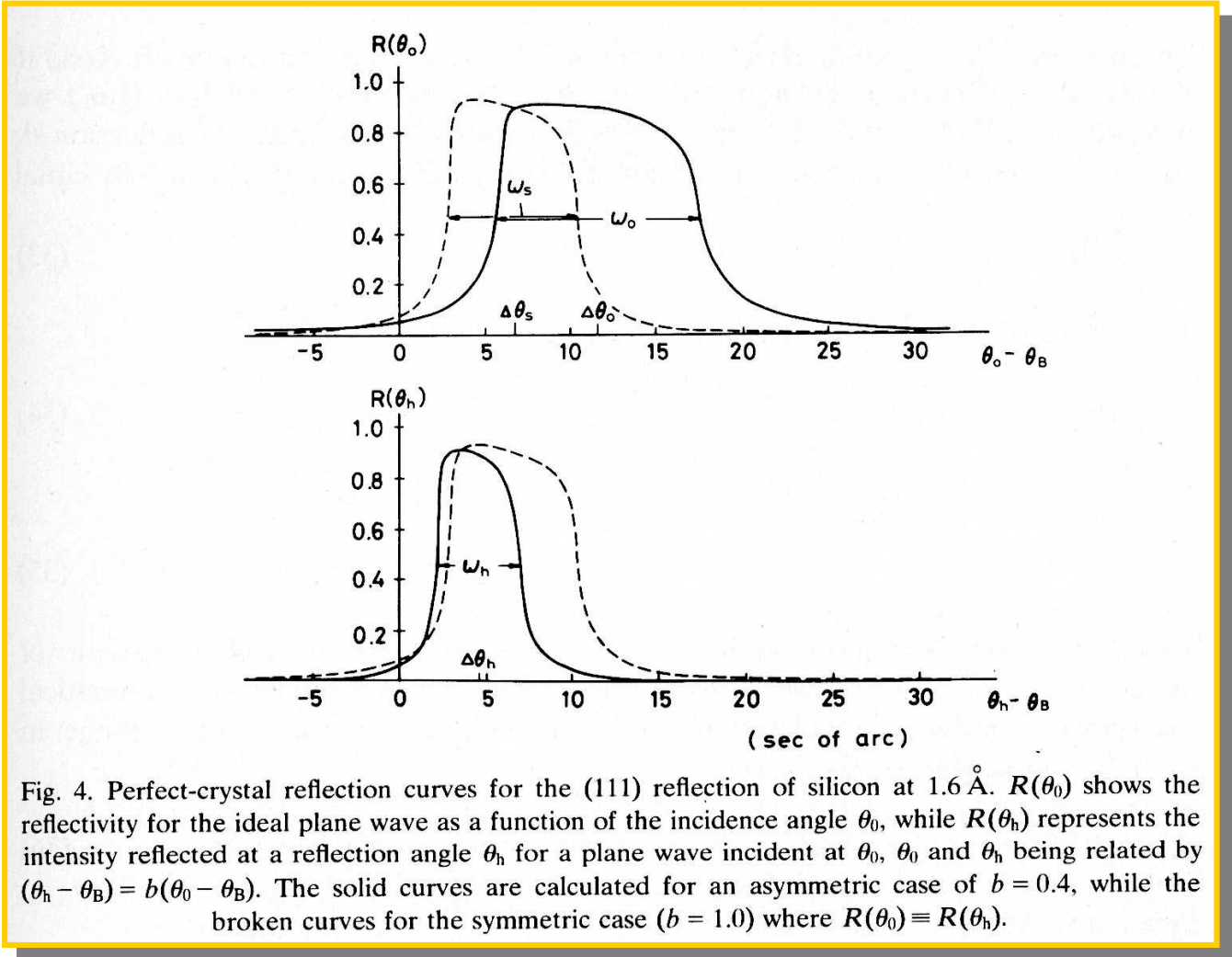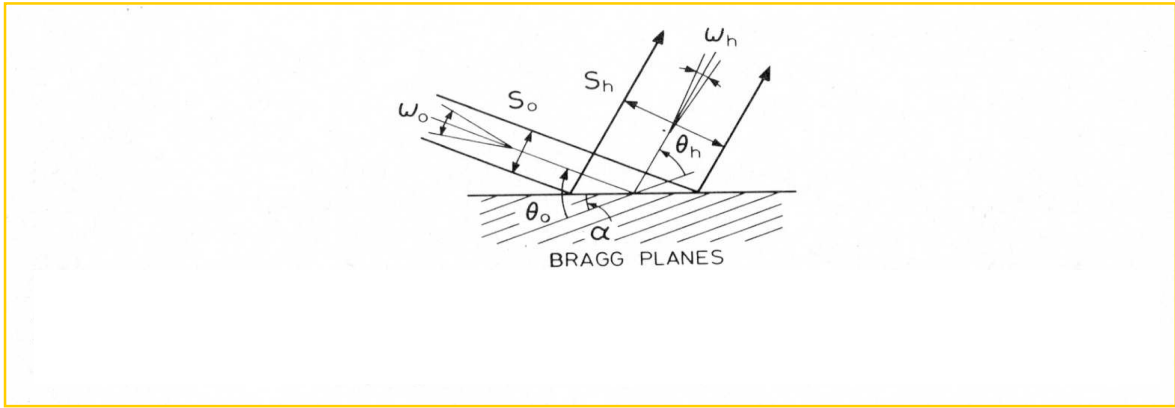

T. Matsushita and H. Hashizume *X-Ray Monochromators*
 Handbook on Synchrotron Radiation, Vol. 1, edited by E.E. Kock
 North-Holland Publishing Company, 1983



Sincrotrone Trieste S.c.p.A.
 Hard-X Ray Optics
 Laboratory


Table 2
Intrinsic Bragg reflection widths ω_s , energy resolutions $\Delta E/E$ and
integral reflecting powers I of perfect crystals of silicon, germanium
and α -quartz at 1.54 Å.

Crystal	hkl	ω_s (second or arc)	$\Delta E/E$ ($\times 10^5$)	I ($\times 10^6$)
Silicon	111	7.395	14.1	39.9
	220	5.459	6.04	29.7
	311	3.192	2.90	16.5
	400	3.603	2.53	19.3
	331	2.336	1.44	11.8
	422	2.925	1.47	15.5
	333	1.989	0.88	9.9
	(511)	1.989	0.88	9.9
	440	2.675	0.96	14.0
531	1.907	0.60	9.3	
Germanium	111	16.338	32.64	85.9
	220	12.449	14.46	67.4
	311	7.230	6.92	37.1
	400	7.951	5.94	42.3
	331	5.076	3.34	25.4
	422	6.178	3.34	32.4
	333	4.127	2.00	20.2
	(511)	4.127	2.00	20.2
	440	5.339	2.14	27.5
531	3.719	1.33	17.7	
α -quartz	100	3.798	10.00	18.8
	101	7.453	15.26	40.9
	110	2.512	3.69	12.2
	10 $\bar{2}$	2.488	3.36	12.9
	200	2.252	2.81	11.5
	112	2.927	3.03	15.5
	202	2.072	1.93	10.6
	212	2.042	1.47	10.7
	20 $\bar{3}$	2.430	1.74	12.9
	301	2.368	1.69	12.6


T. Matsushita and H. Hashizume X-Ray Monochromators
Handbook on Synchrotron Radiation, Vol. 1



Sincrotrone Trieste S.c.p.A.
 Hard-X Ray Optics
 Laboratory

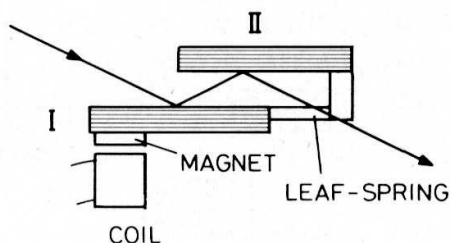

(a)

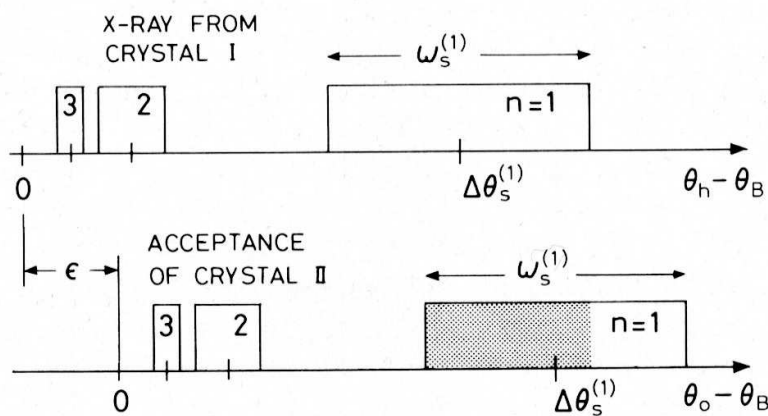

(b)

Fig. 33. An off-set harmonics-rejection monochromator. (a) Geometry of the monochromator. (b) The principle of harmonics rejection. Perfect-crystal reflection curves for the fundamental ($n = 1$) and the harmonics ($n = 2, 3$) are approximated by rectangular boxes. ϵ : off-set or misalignment angle. The shaded area represents delivered X-rays (Hart and Rodrigues 1978).

T. Matsushita and H. Hashizume X-Ray Monochromators Handbook on Synchrotron Radiation, Vol. 1



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory



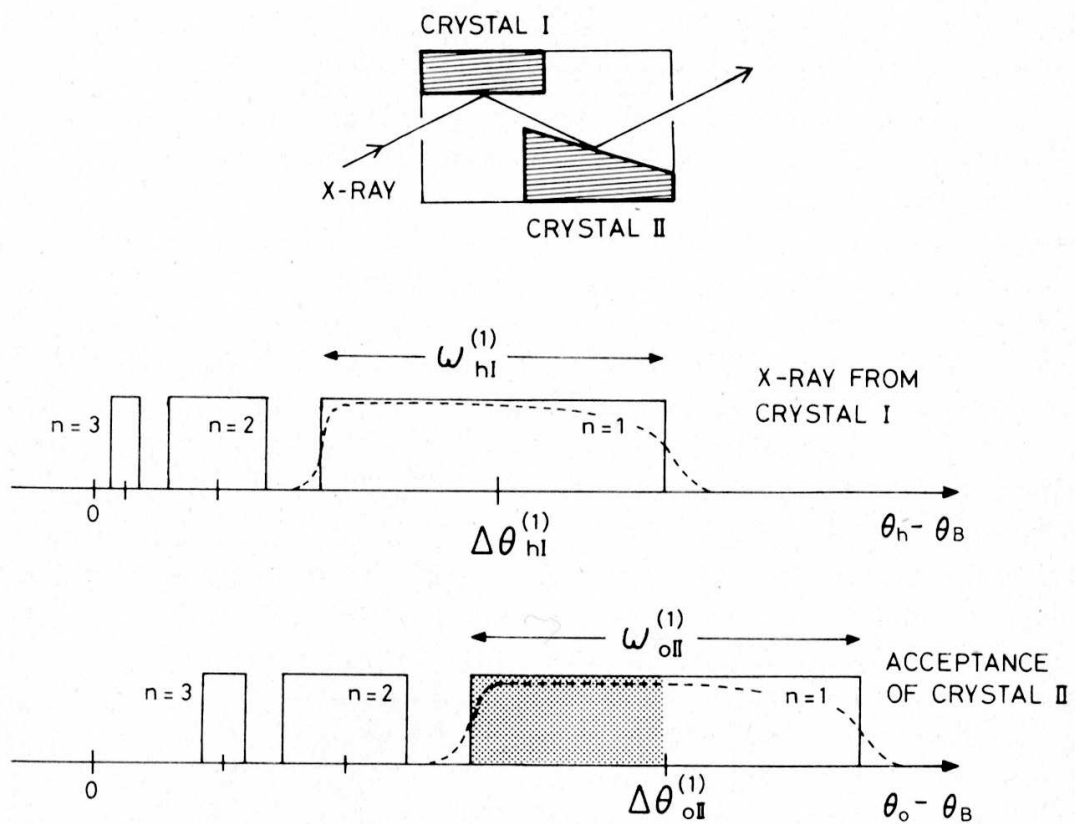


Fig. 34. A monolithic harmonics-rejection monochromator. (a) Crystals I and II of unequal asymmetry factors are built as two outstanding parts of a perfect single crystal. (b) The principle of harmonics rejection. Perfect-crystal reflection curves for the fundamental ($n = 1$) and the harmonics ($n = 2, 3$) are approximated by rectangular boxes. The broken curves show the real reflection curves for the fundamental. The shaded area represents delivered X-rays.

T. Matsushita and H. Hashizume *X-Ray Monochromators Handbook on Synchrotron Radiation, Vol. 1*

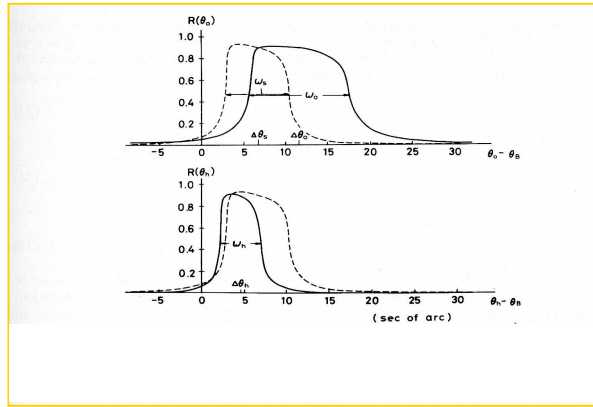

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

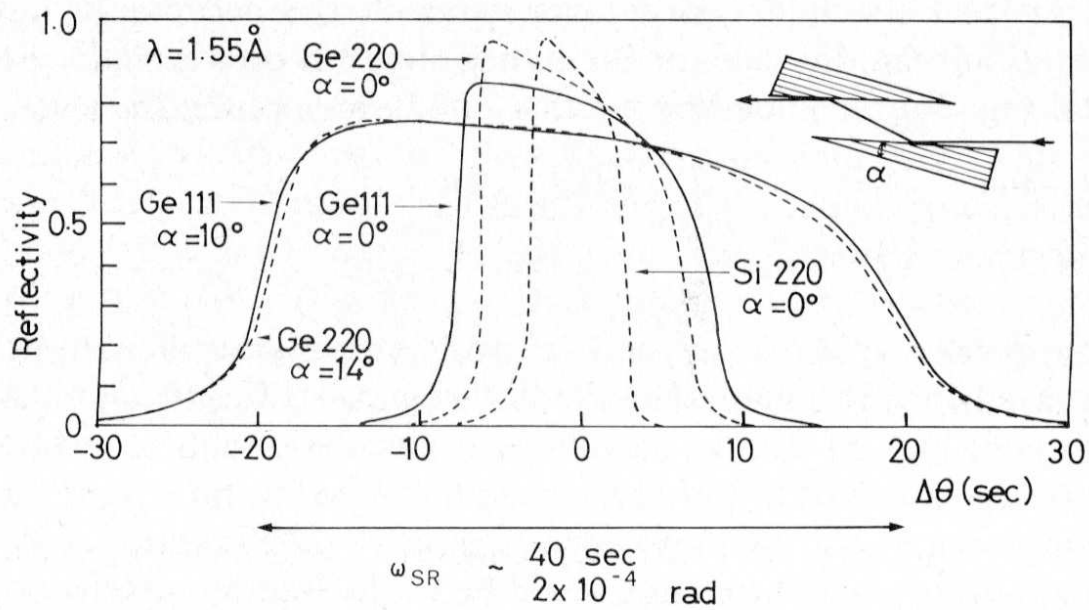

Calculated reflectivity curves of grooved monochromators using various asymmetric reflections of silicon and germanium for 1.55 Å X-rays (Kohra et al. 1978).


T. Matsushita and H. Hashizume X-Ray Monochromators Handbook on Synchrotron Radiation, Vol. 1

Note as the refractive effect on the first crystal has been totally compensated by the second one





Sincrotrone Trieste S.c.p.A.
 Hard-X Ray Optics
 Laboratory

Double crystal monochromator

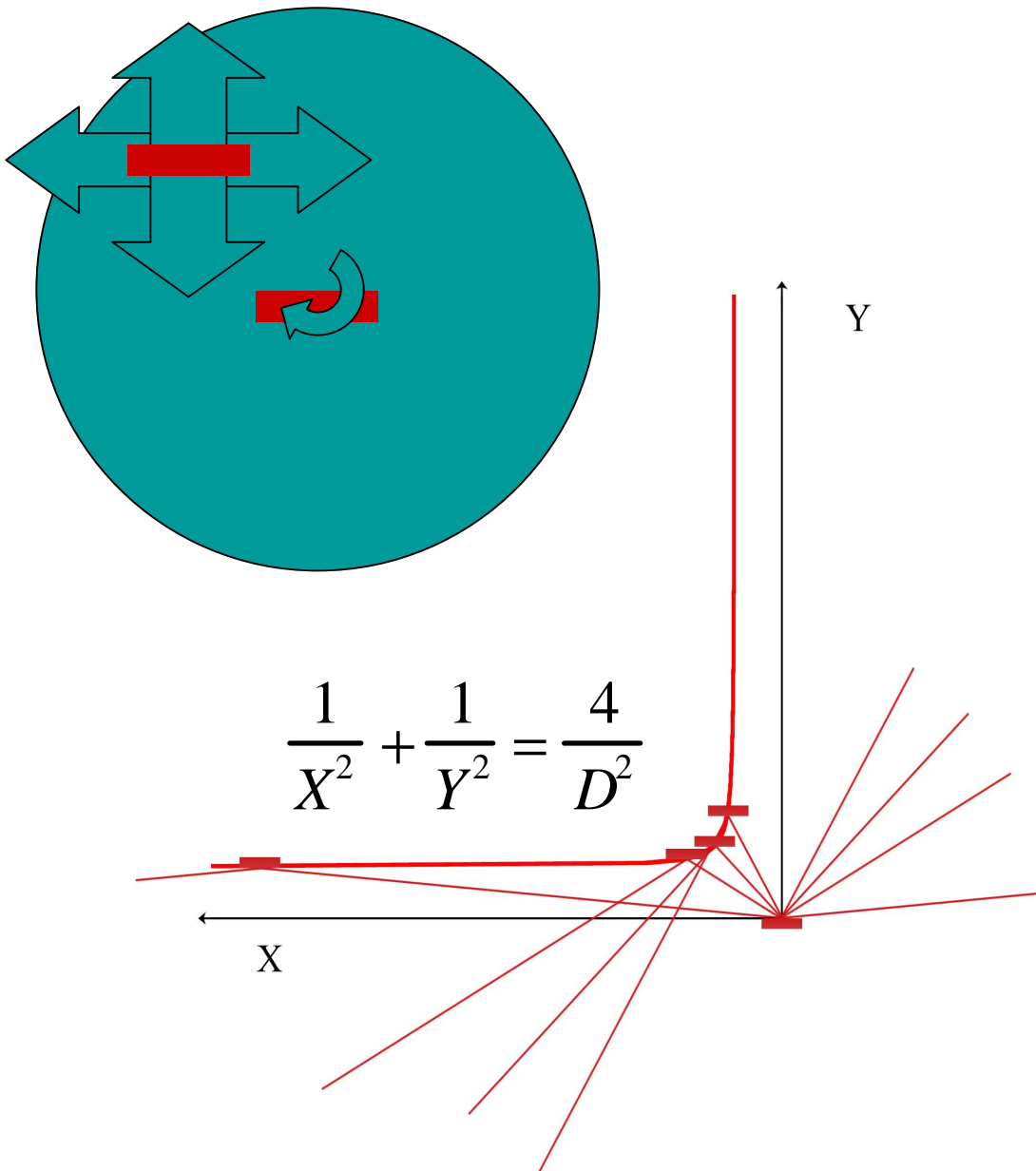


Double crystal monochromator

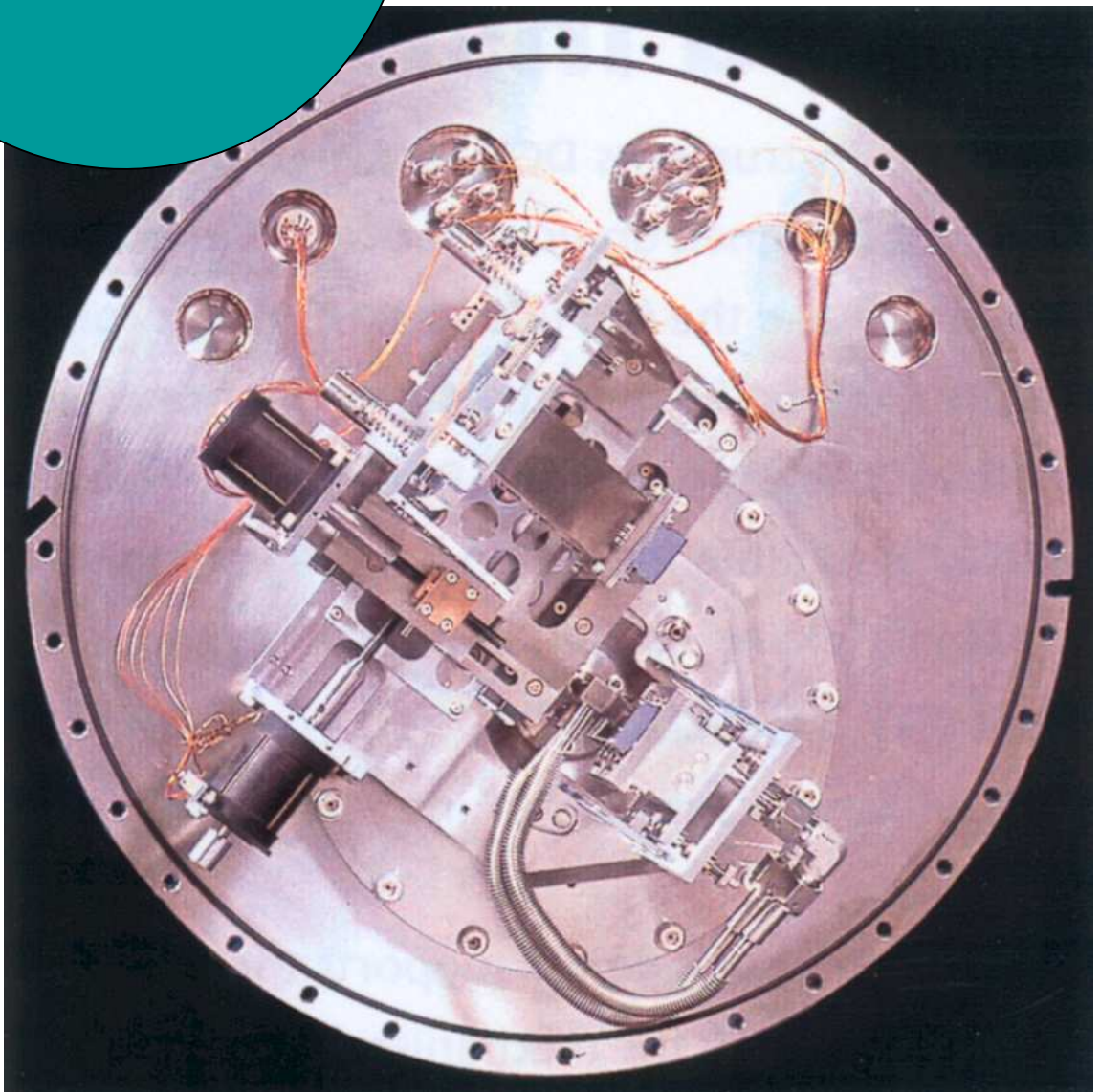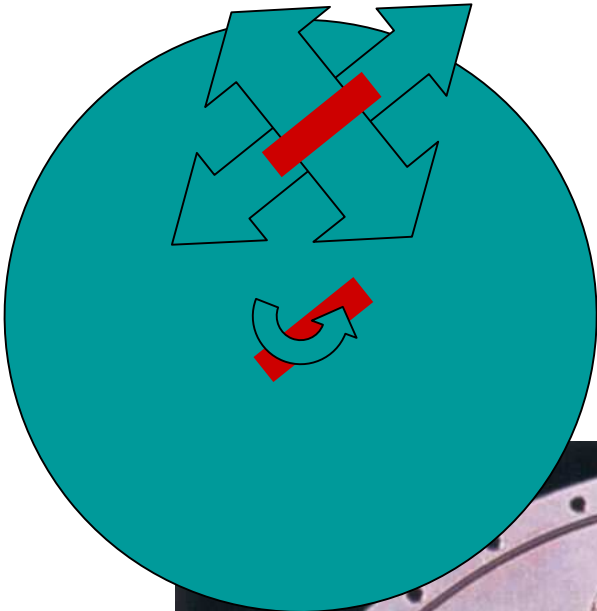

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory


Diffraction 1

57 poles wiggler source at ELETTRA
400mA, 1.6T and 2GeV
total power: 8 kW

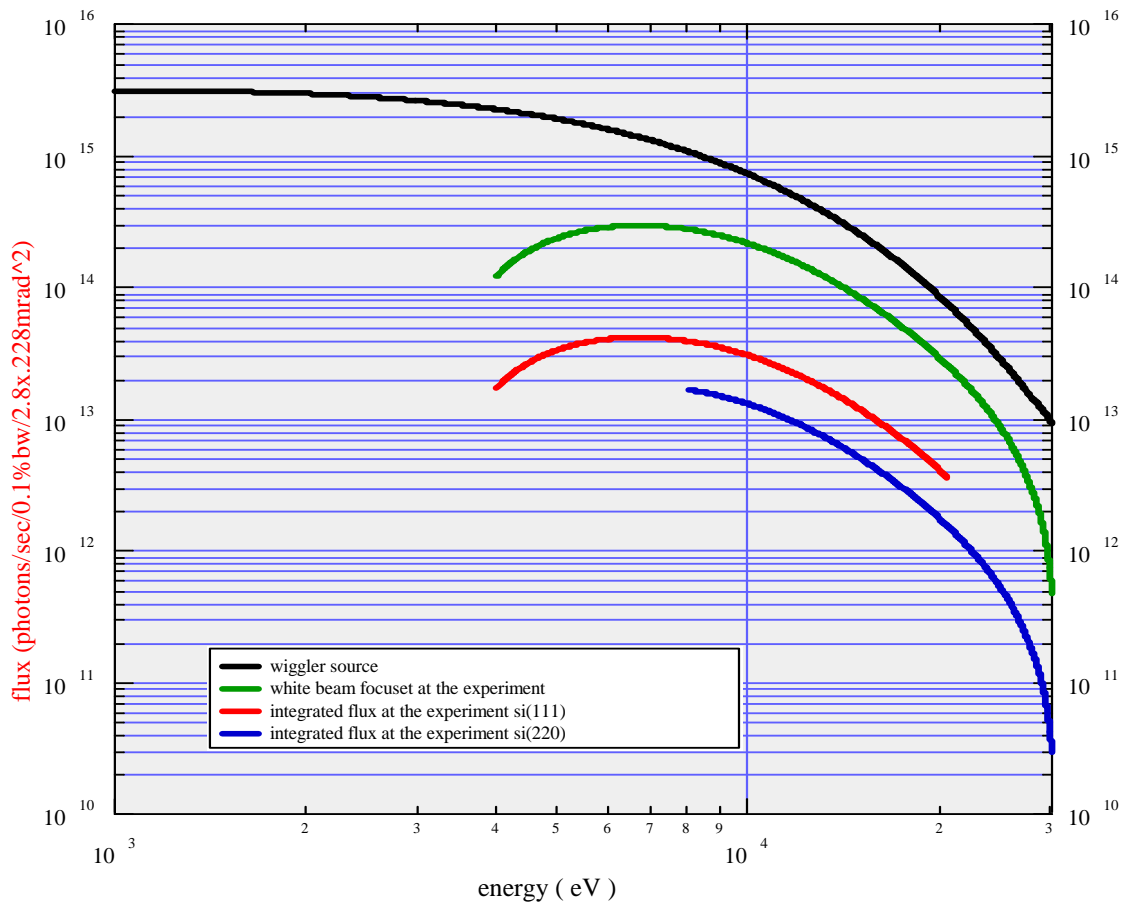

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

The double crystal
monochromator at the **diffraction1** beamline

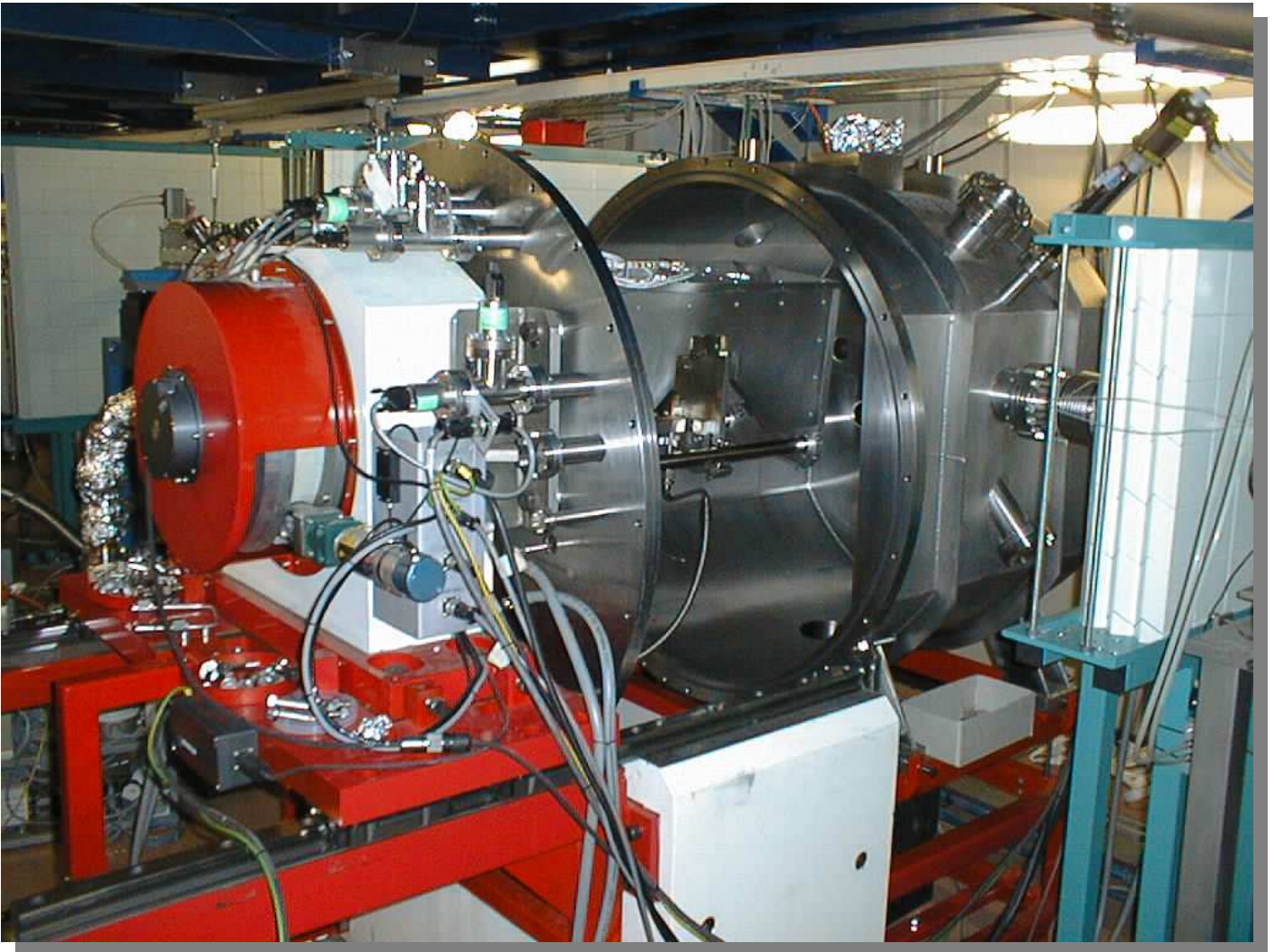

*Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory*

Diffraction 1:

first optical element in the beam

Si(111) internally water cooled


total power absorbed 0.5 kW*

* $1.5 \times 0.28 \text{ mrad}^2$




*Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory*

Conceptual design of an internally water cooled crystal

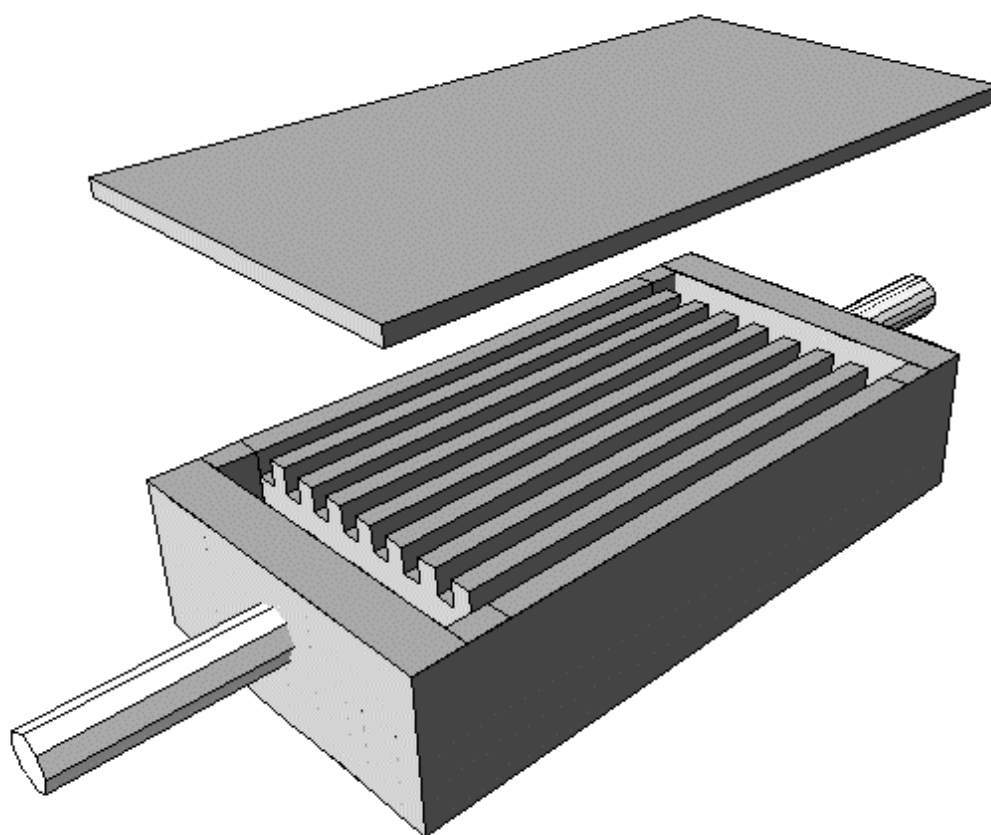

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

The two Si components before the Si-Si brazing

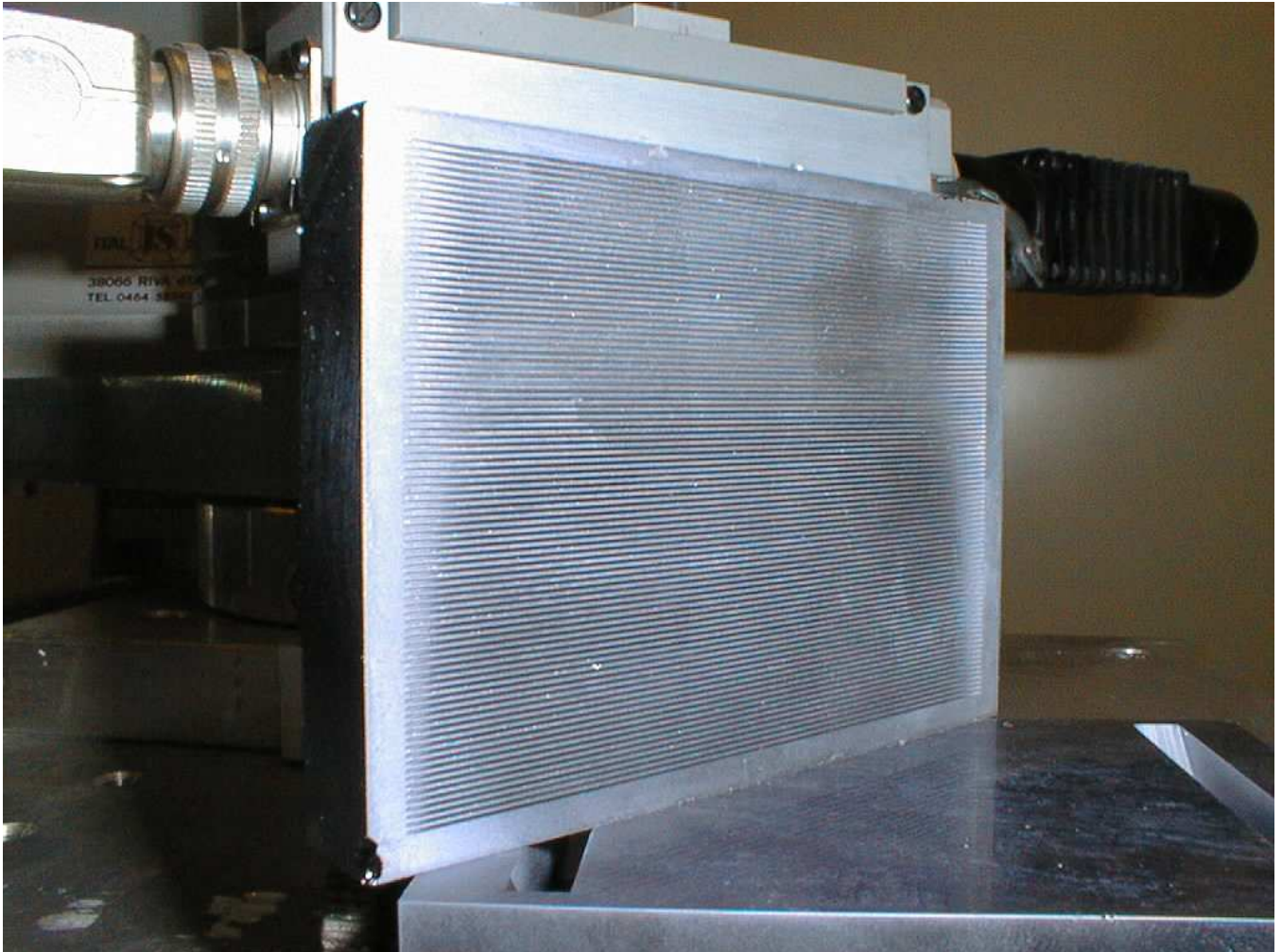

channels:


thickness: 300 μ m

depth: 2mm



Elettra
Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

3D picture of the double ϑ - 2ϑ equipment

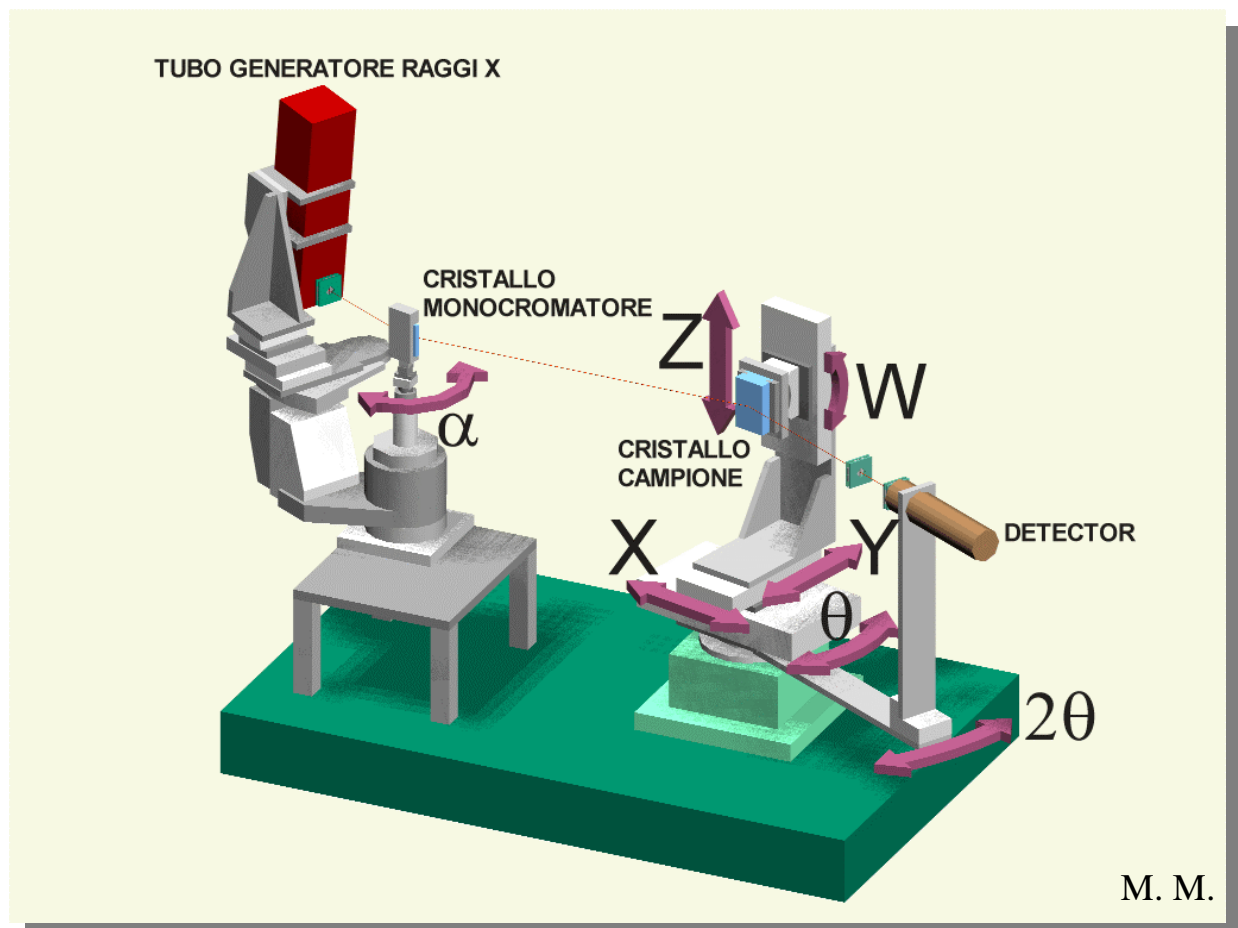

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

The double ϑ - 2ϑ test station at the Hard X-ray Laboratory



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory


Scanning panel



Scanning parameters

Current Axis	Current position	
THETA	:	-47.4802 [deg]
Center	Amplitude	Steps
-47.4820 [deg]	0.0100 [deg]	51
Counting time:	2 s	SET TIME...

Scanning graph for current axis

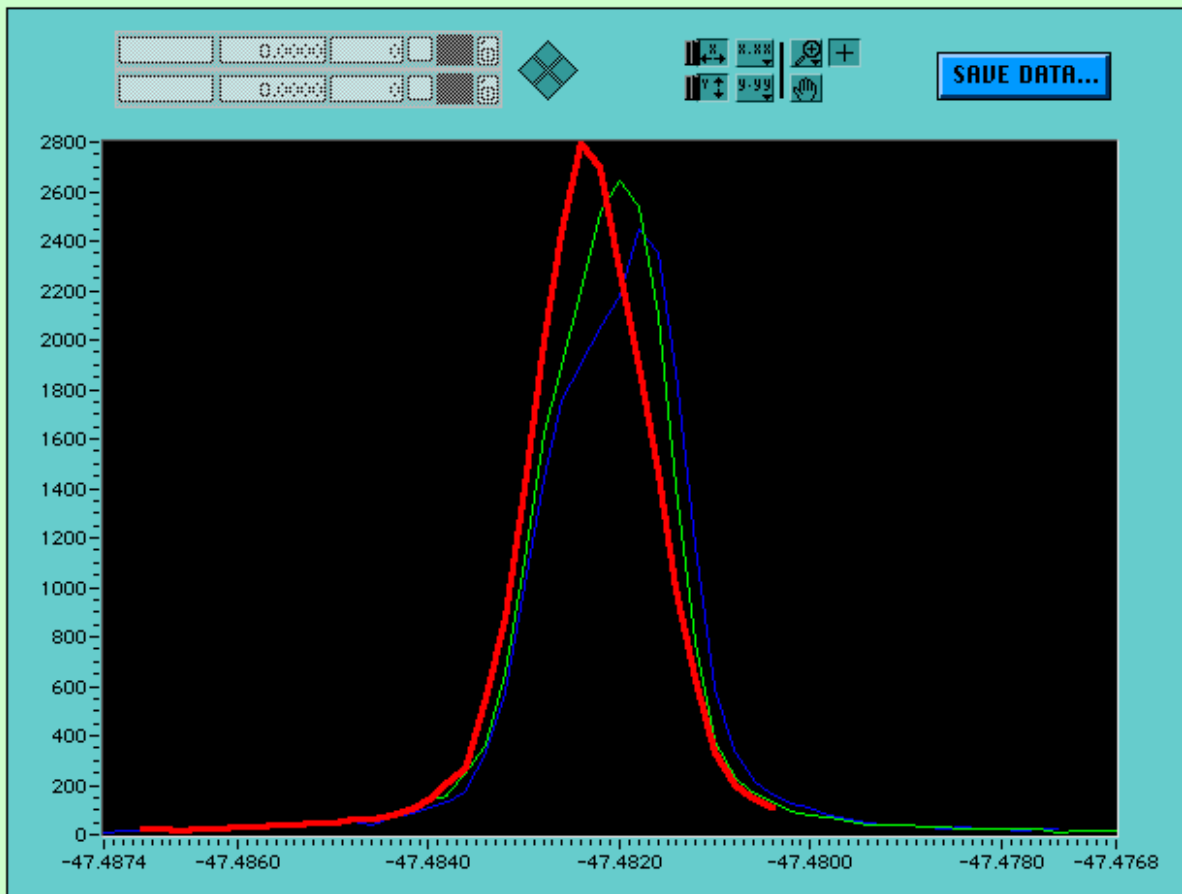


Action commands

SCAN
STOP
MOVE TO...
Exit

Scanning on axis SMC9000.1 (step 35 of 51)



DISTRIBUZIONE DELL'AMPIEZZA DELLE ROCKING-CURVE VALUTATE MEDIANTE DIFFRAZIONE SUI PIANI Si-111

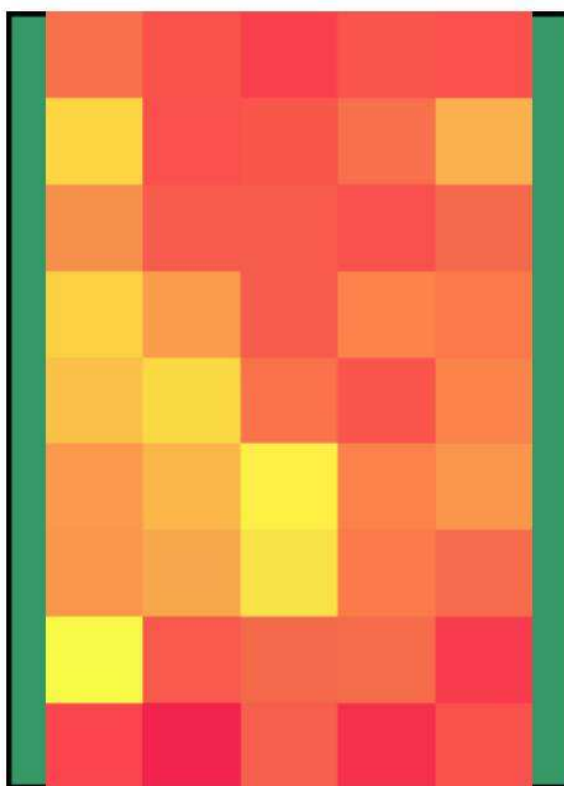


14.1
13.6
13.1
12.6
12.1
11.6
11.1
10.7
10.2
9.7

Valori espressi in
secondi di grado




Topography of the internal cooled Si-crystal with channels perpendicular to the scattering plane

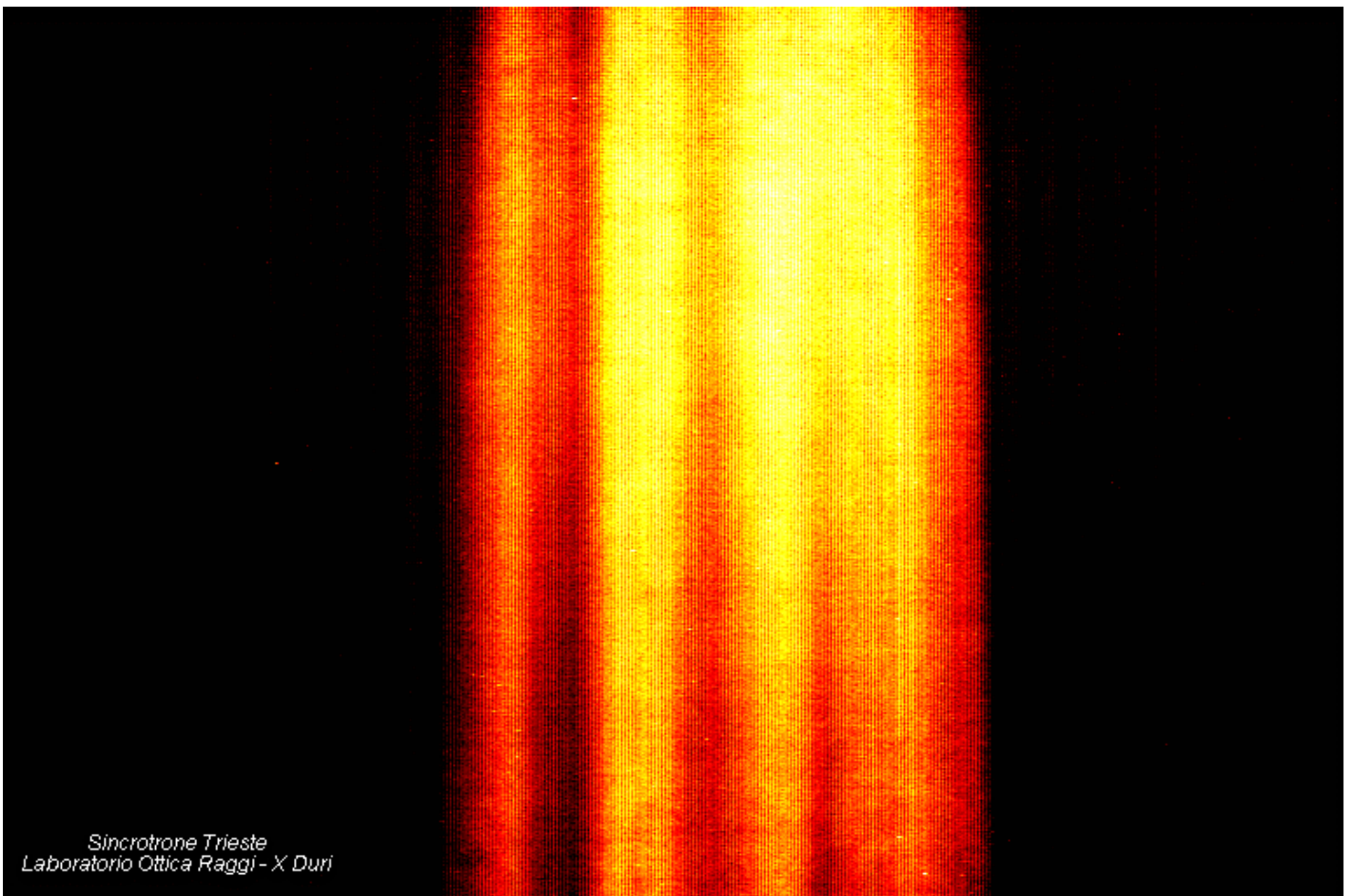


*Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory*

Same topography but with
channels in the same direction of
the scattering plane

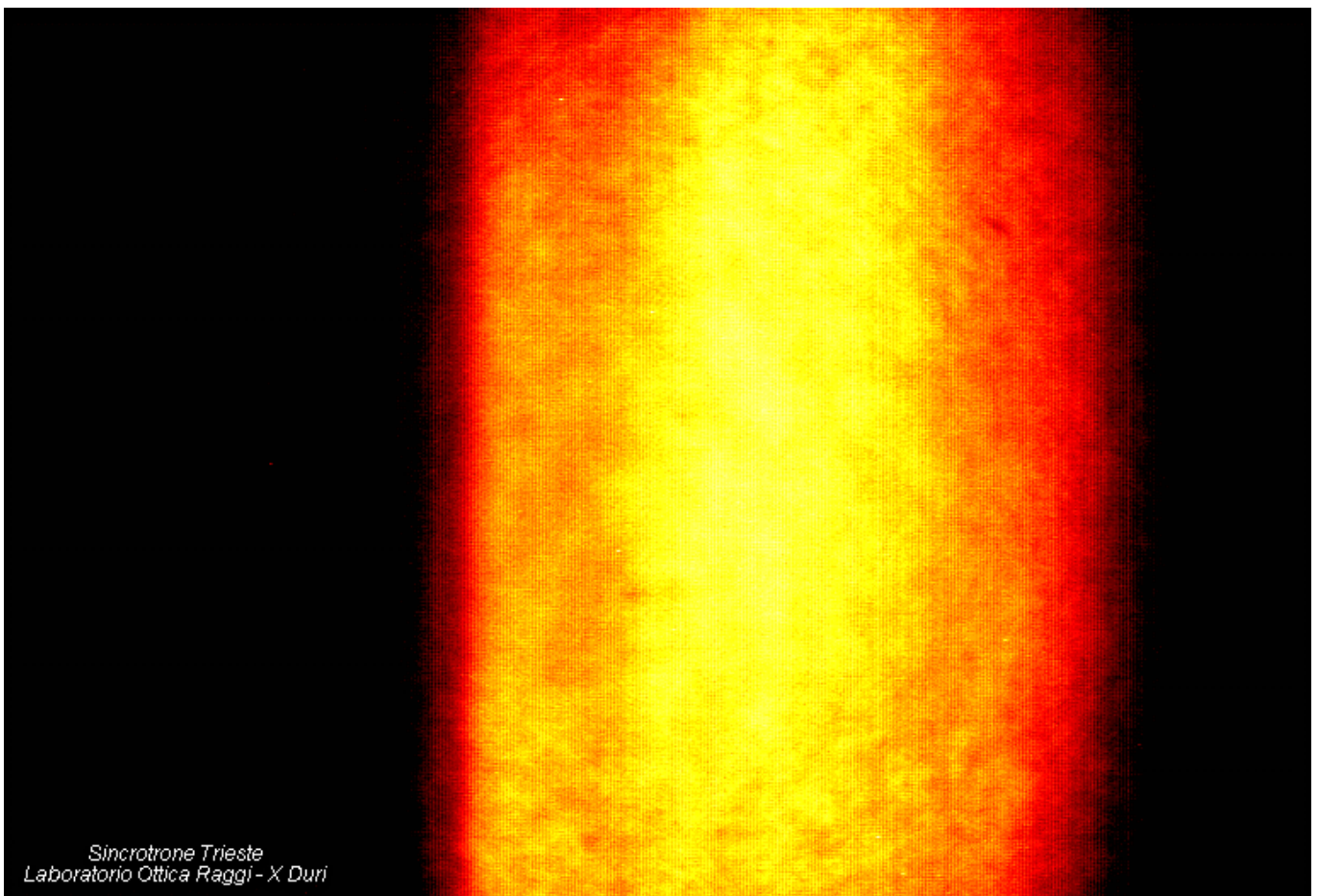


Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

surface of a Si(111) crystal with a evident stressed structure induced by a back-side machining and not removed by chemical hatching

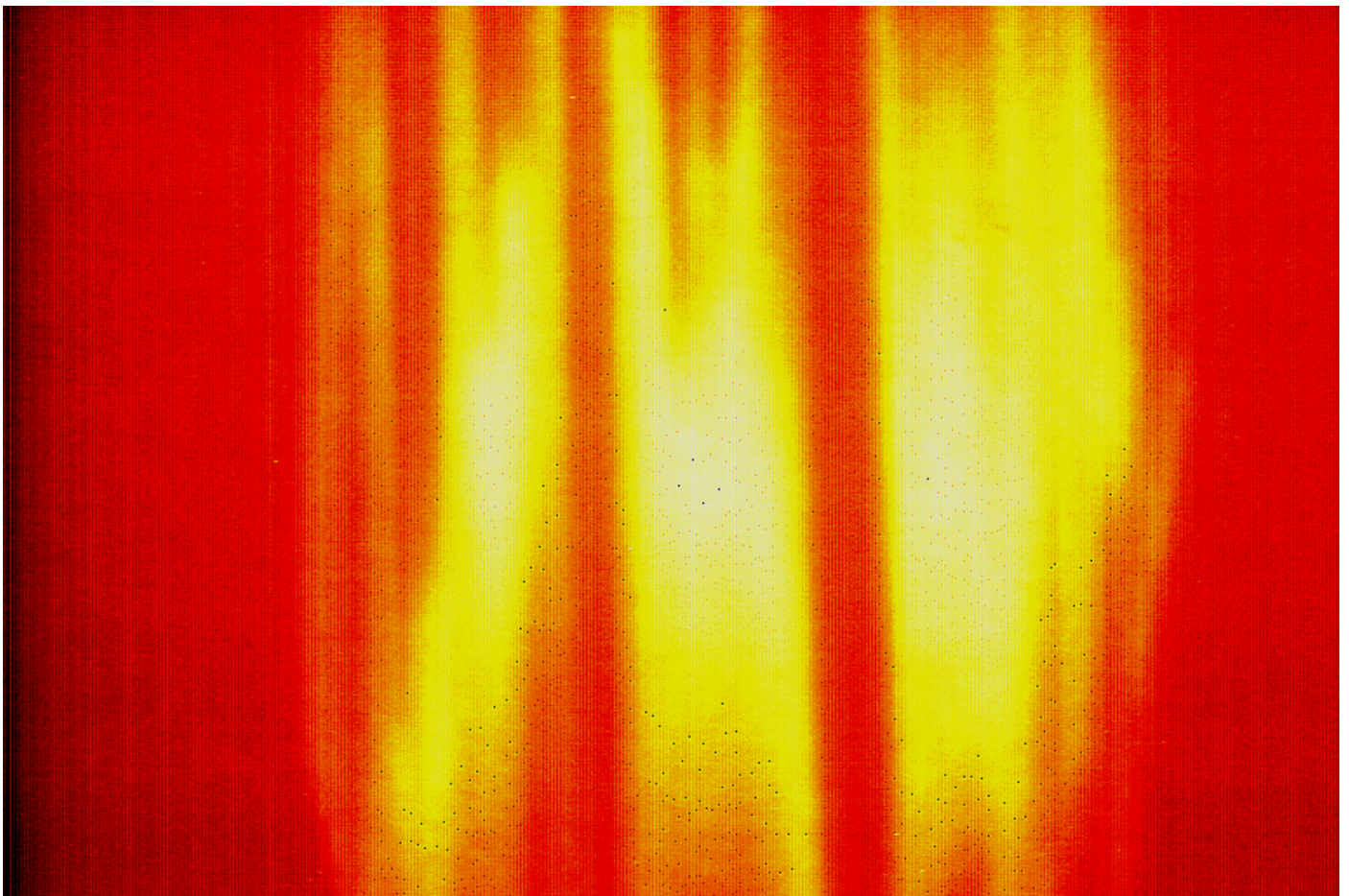



*Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory*

The inclined double crystal monochromator setup to reduce the power density

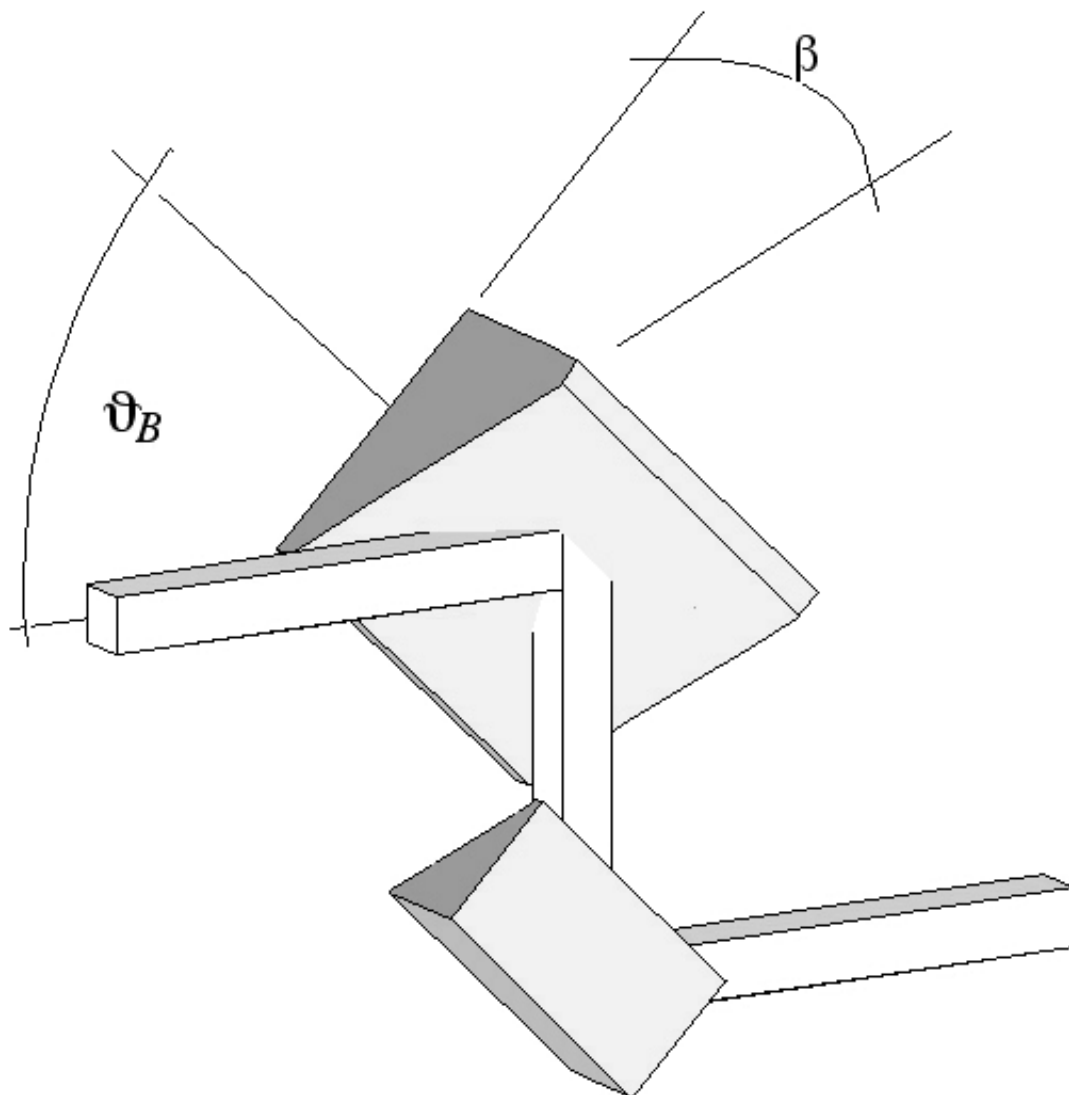

Si(111) inclined channel-cut crystal monochromator designed for the ALOISA beamline.

energy range: 2.8 to 8 KeV
beam dimension: 3x3 mm
source: wiggler-ondulator

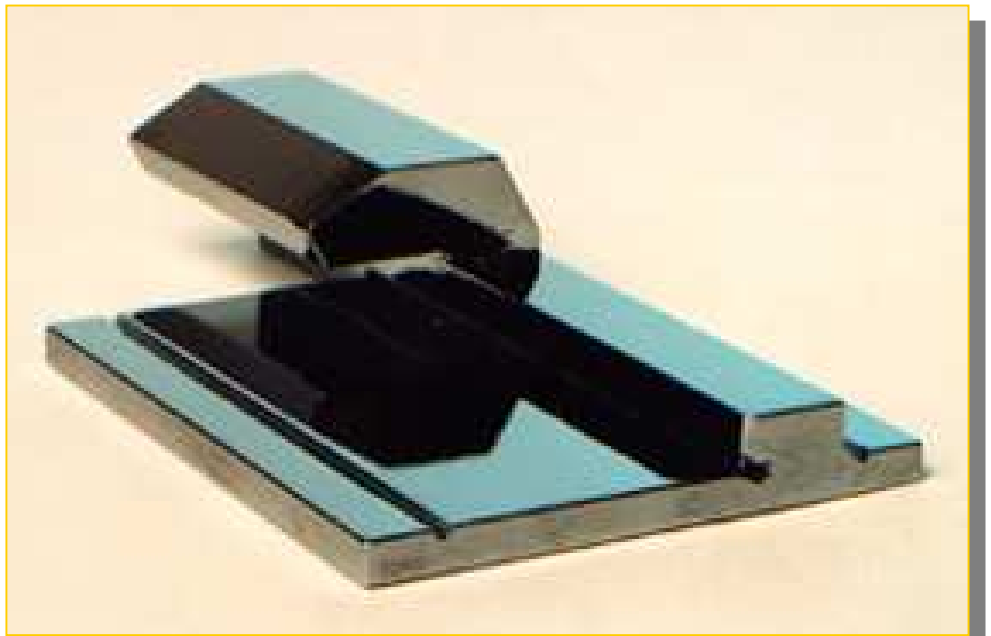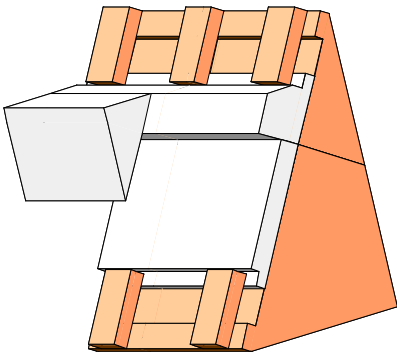


E. Busetto et al.: *"The High Energy Monochromator for the ALOISA Beamline"*. Rev. Sci. Instrum. **66** (2), February 1995



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

*A new prototype of monochromator
under test with x-rays*

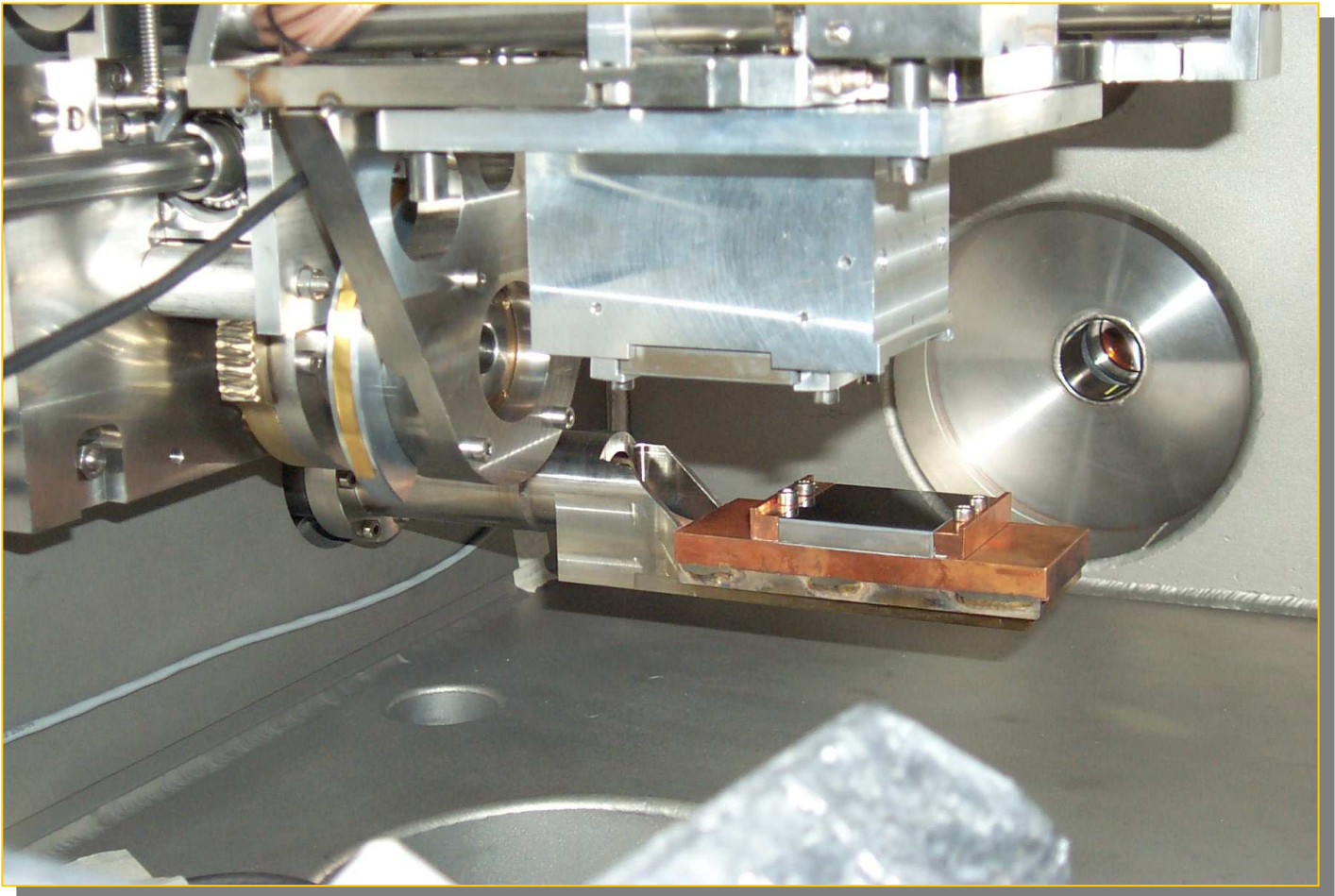

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

- Detector for Hard X-rays

- Two large families:

single counters

integrators




- Single counters :
- These systems allow to collect all the electrons produced by the absorption of an x-ray.
- The mean number of electrons produced during the absorption process is proportional to energy of the single x-ray.



Scintillatori

I rivelatori a scintillazione sono il risultato dell'accoppiamento di un cristallo isolante drogato e di un convertitore-moltiplicatore optoelettronico, il fototubo

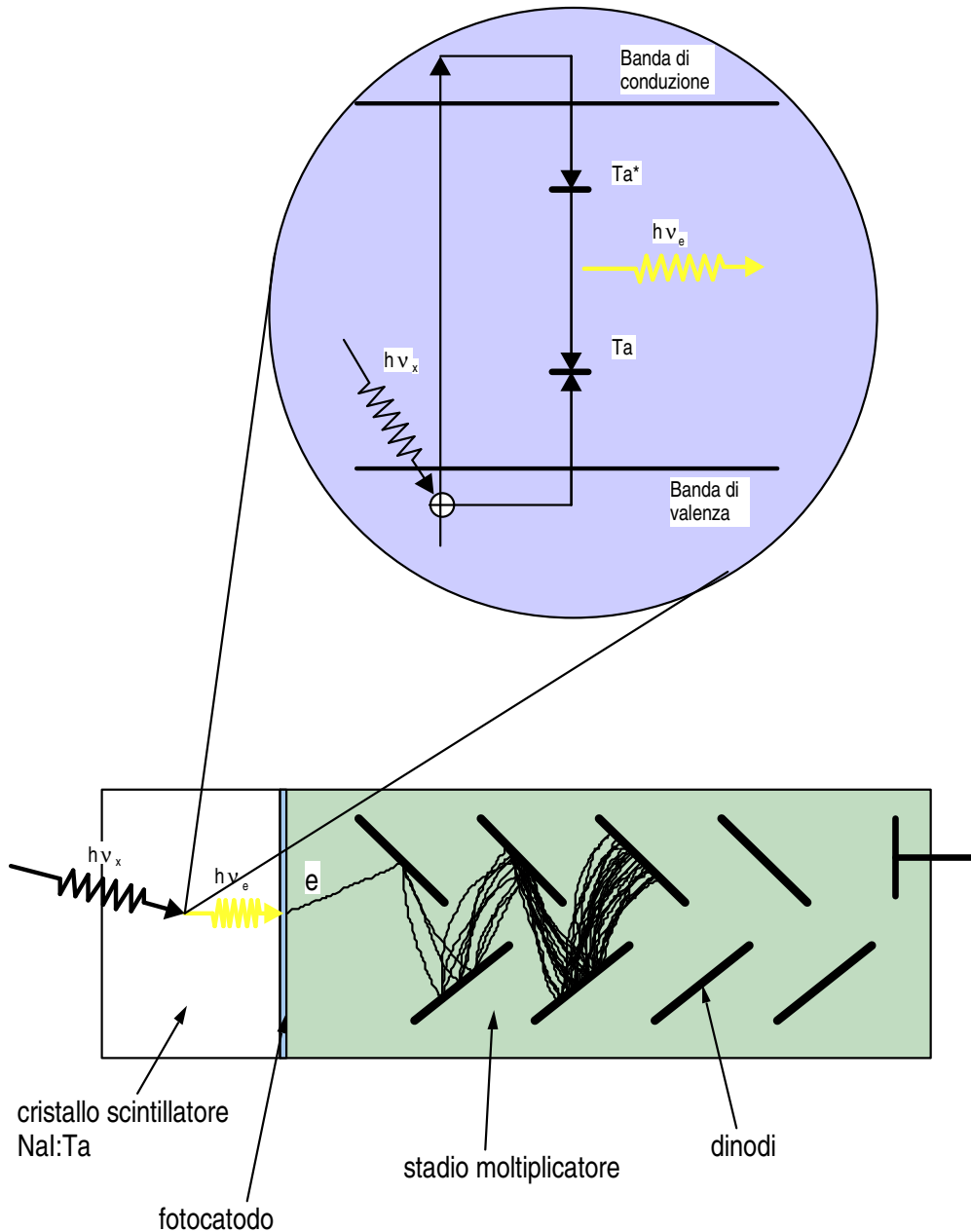

Integratori

- Sono sistemi di rivelazione “integrata” nel tempo dove si perde la correlazione diretta fra elettroni prodotti ed energia del fotone incidente.
- L’intensità del segnale locale deve essere proporzionale al numero di fotoni assorbiti nella stessa zona.



La pellicola radiografica


è stato ed è tuttora uno tra i più usati rivelatori ad integrazione sfrutta la sensibilità alle radiazioni elettromagnetiche degli alogenuri di argento in particolare dell' AgBr. La reazione fotochimica produce ioni di Ag^+ con una densità proporzionale alla quantità di radiazione assorbita. Le sostanze fortemente riducenti utilizzate negli sviluppi trasformano Ag^+ , localmente prodotto, in Ag metallico dalla tipica colorazione nera.




Elettra
Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Curva caratteristica della densità contro il tempo di esposizione

Linearità e range dinamico delle emulsioni fotografiche

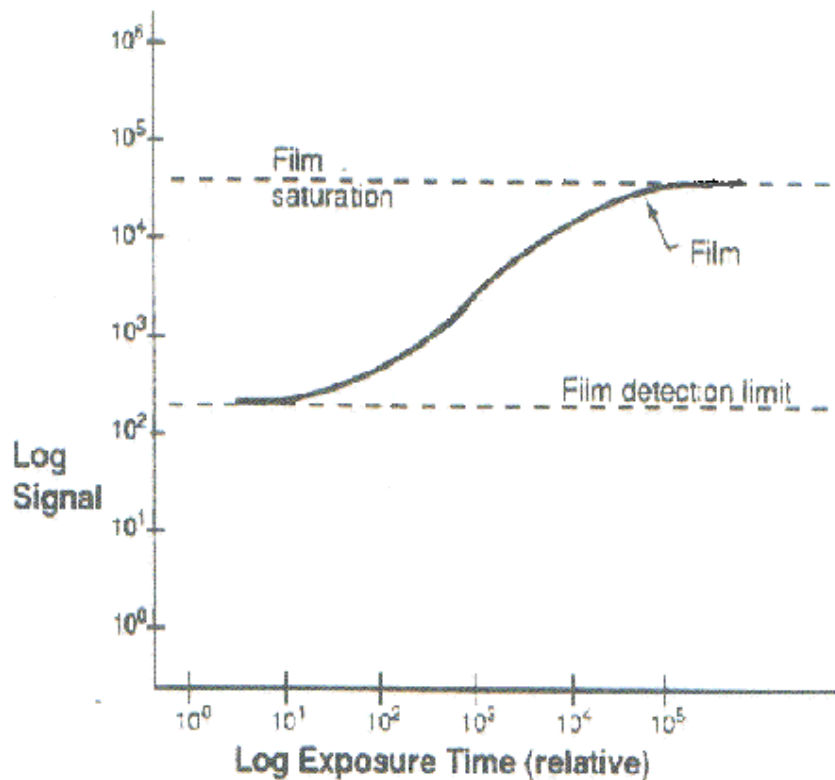

Integratori bidimensionali digitali

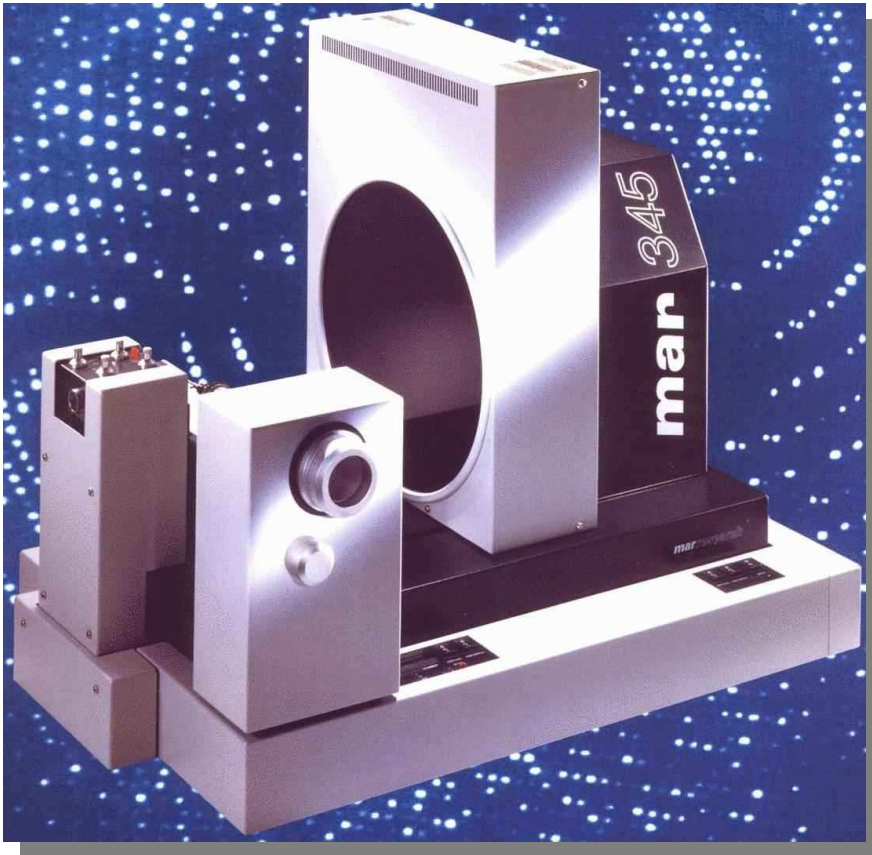


Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Spaccato di un rivelatore CCD commerciale per raggi X





Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Linearità, range dinamico ed efficienza a confronto

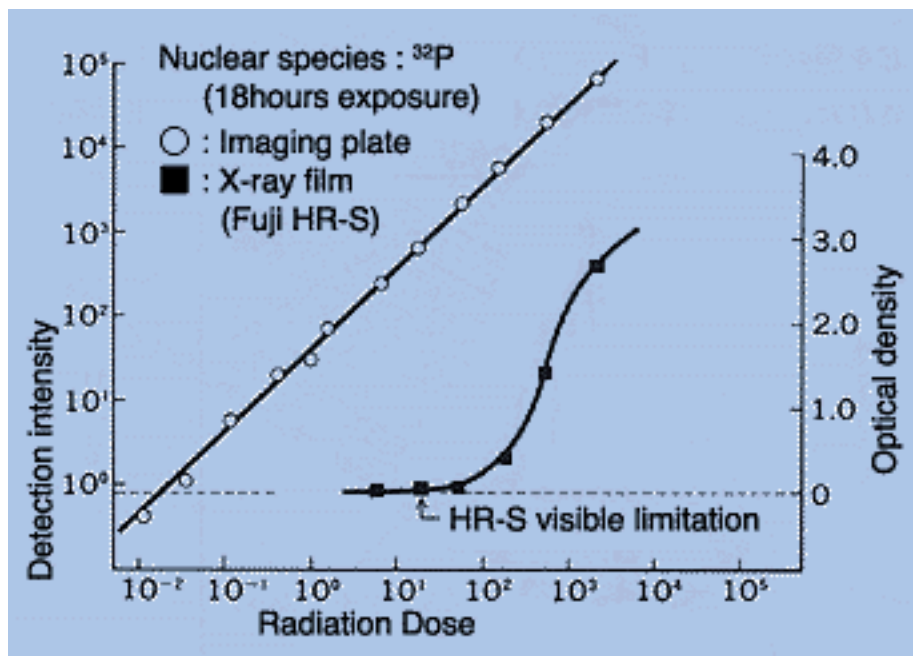


Fisica dell'Image Plate.....un rivelatore analogico a lettura digitale

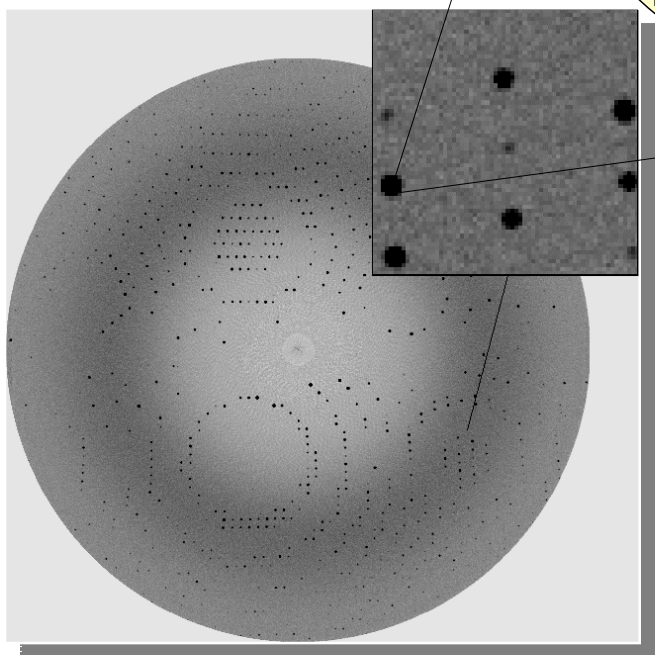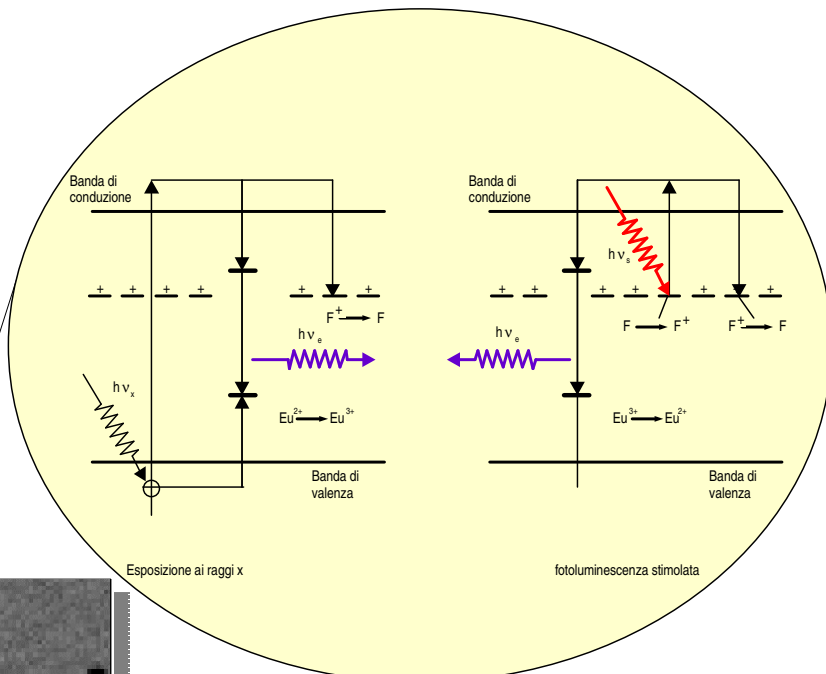

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

I rivelatori a gas : principio di funzionamento

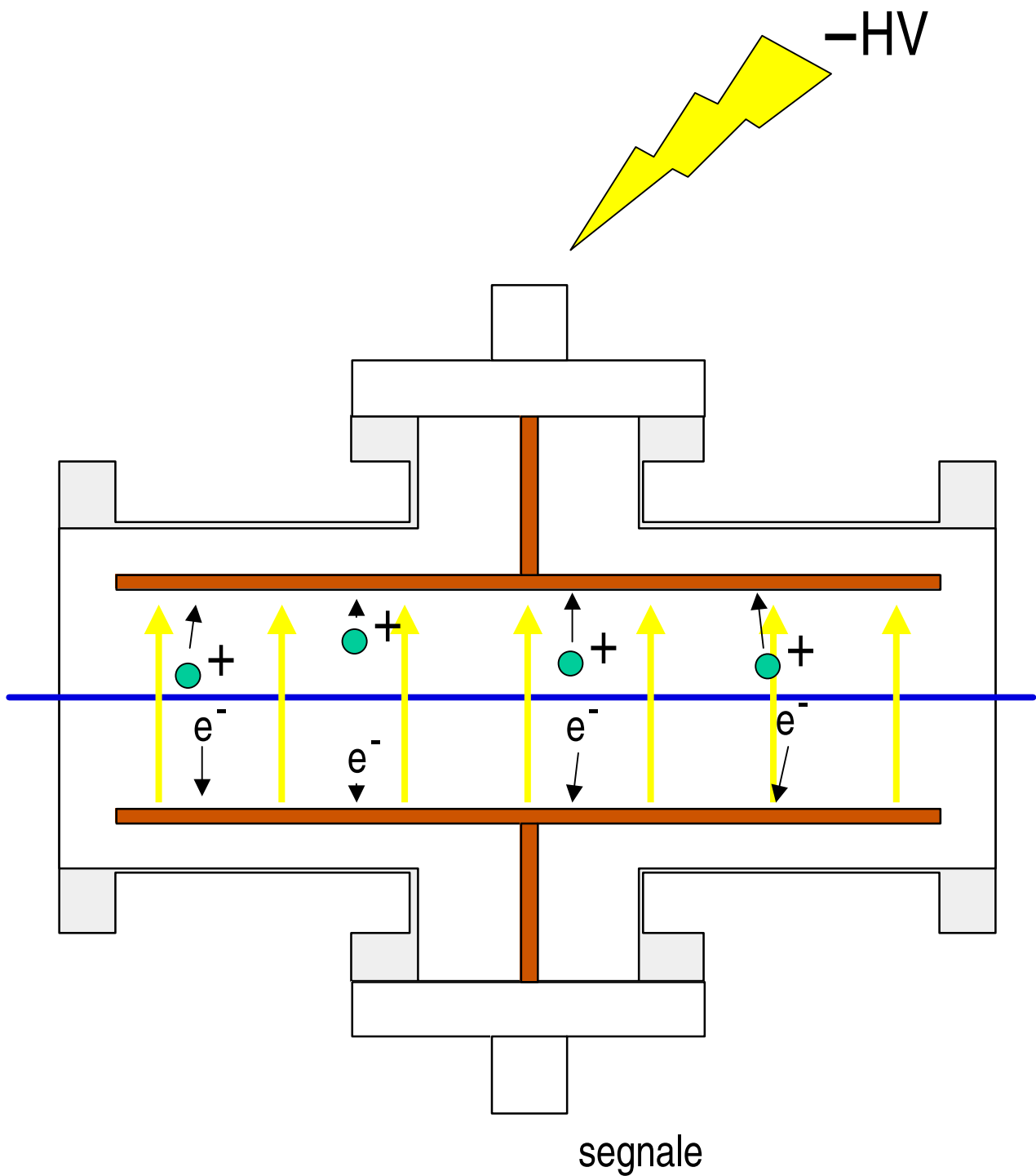

Wiggler insertion devices at ELETTRA

Permanent magnet Wiggler W14.0 XRD1 source:

*	B_0	= 1.6 Tesla
*	Period length	= 140 mm
*	No. of poles	= 59
*	Total length	= 4500 mm
*	K_y	= 19.6
*	E_{cm}	= 4.2 KeV (2GeV) / 6.0KeV (2.4GeV)


Multipole Superconducting Wiggler XRD2 source:

*	B_0	= 3.5 Tesla
*	Period length	= 64 mm
*	No. of poles	= 49
*	Total length	= 1568 mm
*	K_y	= 20.9
*	E_{cm}	= 9.2 KeV(2GeV) / 13.2KeV(2.4GeV)




Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Diffraction2 conceptual layout

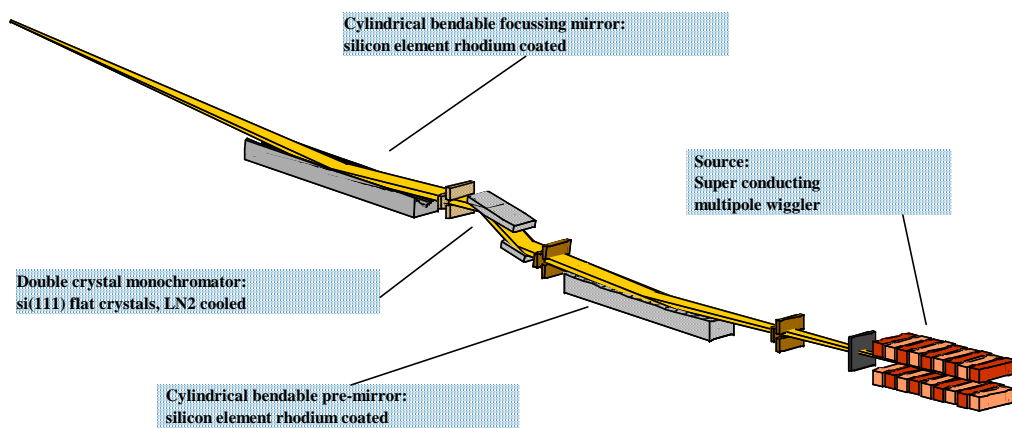

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Diffraction2 conceptual layout

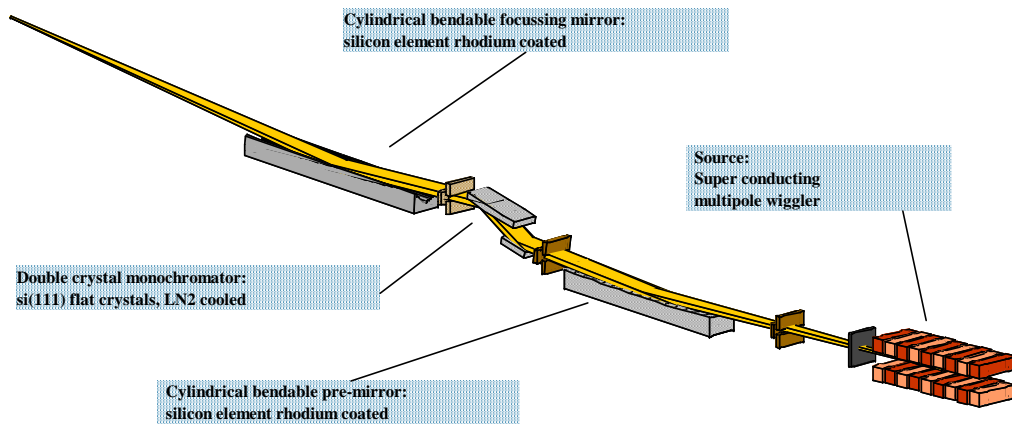

Premirror:

Final shape:

Tangential radius:

Grazing angle:

Source distance:

cylindrical vertical collimator

13.5 Km

0.18° (200 μ rad vertical)

21139 mm


Active optical surface:

1400 x 45 mm²



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Diffraction2 conceptual layout

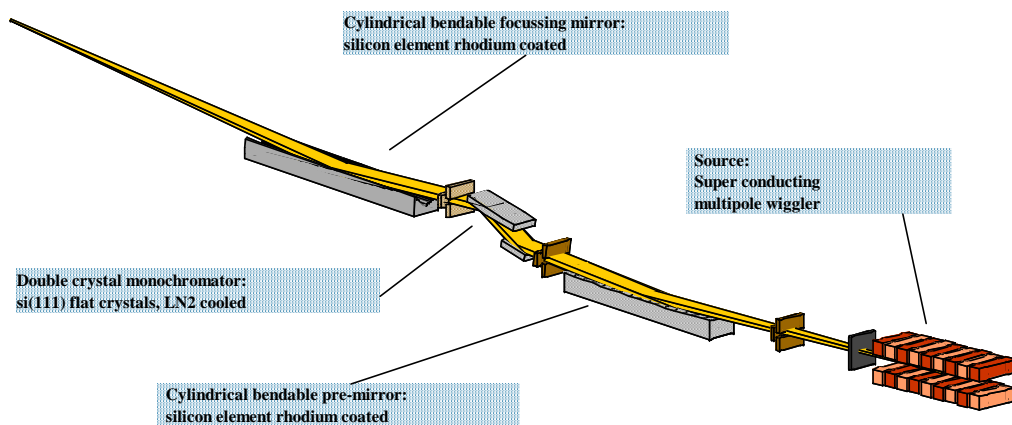

Premirror:

Final shape:

Tangential radius:

Grazing angle:

Source distance:

cylindrical vertical collimator

13.5 Km

0.18° (200 μ rad vertical)

21139 mm

Active optical surface:

1400 x 45 mm²

Focussing mirror:

Final shape

toroidal

Tangential radius:

6.0 Km (0.18°)

Sagittal radius:

49.9 mm (0.18°)

Grazing angle:

0.18°

Source distance:

26932 mm

Active optical surface:

1400 x 55 mm²

Focal distance:

11000 mm with 0.18°

Demagnification:


1.9:1 vertical plane

2.4:1 horizontal plane



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

Diffraction2 conceptual layout

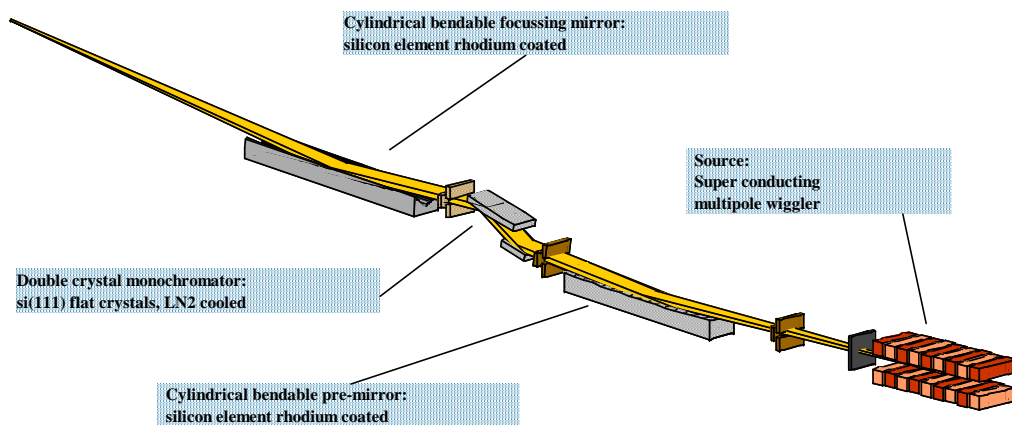

Premirror:

Final shape:

Tangential radius:

Grazing angle:

Source distance:

cylindrical vertical collimator

13.5 Km

0.18° (200 μ rad vertical)

21139 mm

Active optical surface:

1400 x 45 mm²

Focussing mirror:

Final shape

toroidal

Tangential radius:

6.0 Km (0.18°)

Sagittal radius:

49.9 mm (0.18°)

Grazing angle:

0.18°

Source distance:

26932 mm

Active optical surface:

1400 x 55 mm²

Focal distance:

11000 mm with 0.18°

Demagnification:

1.9:1 vertical plane

2.4:1 horizontal plane

Horizontal acceptance:

1 mrad max


Vertical acceptance:

200 μ rad



*Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory*

Diffraction2 conceptual layout

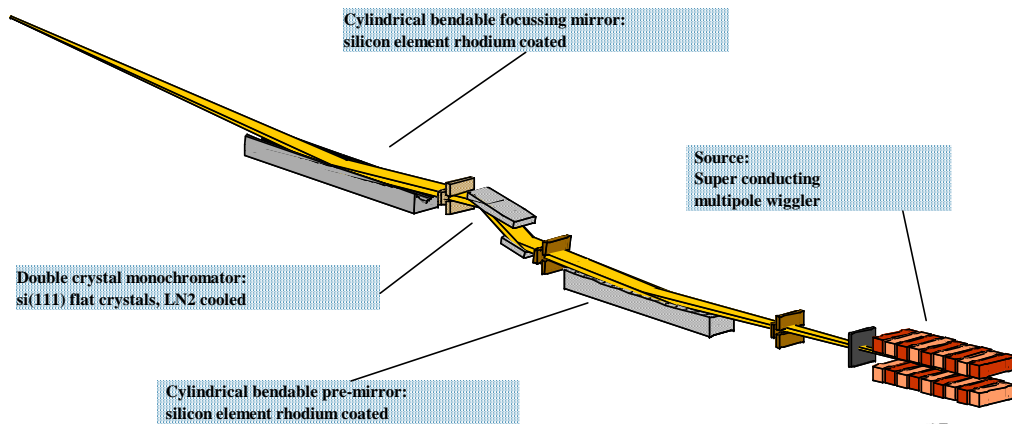

Premirror:


Final shape: cylindrical vertical collimator
 Tangential radius: 13.5 Km
 Grazing angle: 0.18° (200 μrad vertical)
 Source distance: 21139 mm

Active optical surface: 1400 x 45 mm²

Focussing mirror:

Final shape: toroidal
 Tangential radius: 6.0 Km (0.18°)
 Sagittal radius: 49.9 mm (0.18°)
 Grazing angle: 0.18°
 Source distance: 26932 mm
 Active optical surface: 1400 x 55 mm²
 Focal distance: 11000 mm with 0.18°
 Demagnification: 1.9:1 vertical plane
 2.4:1 horizontal plane

Horizontal acceptance: 1mrad max
Vertical acceptance: 200 μrad

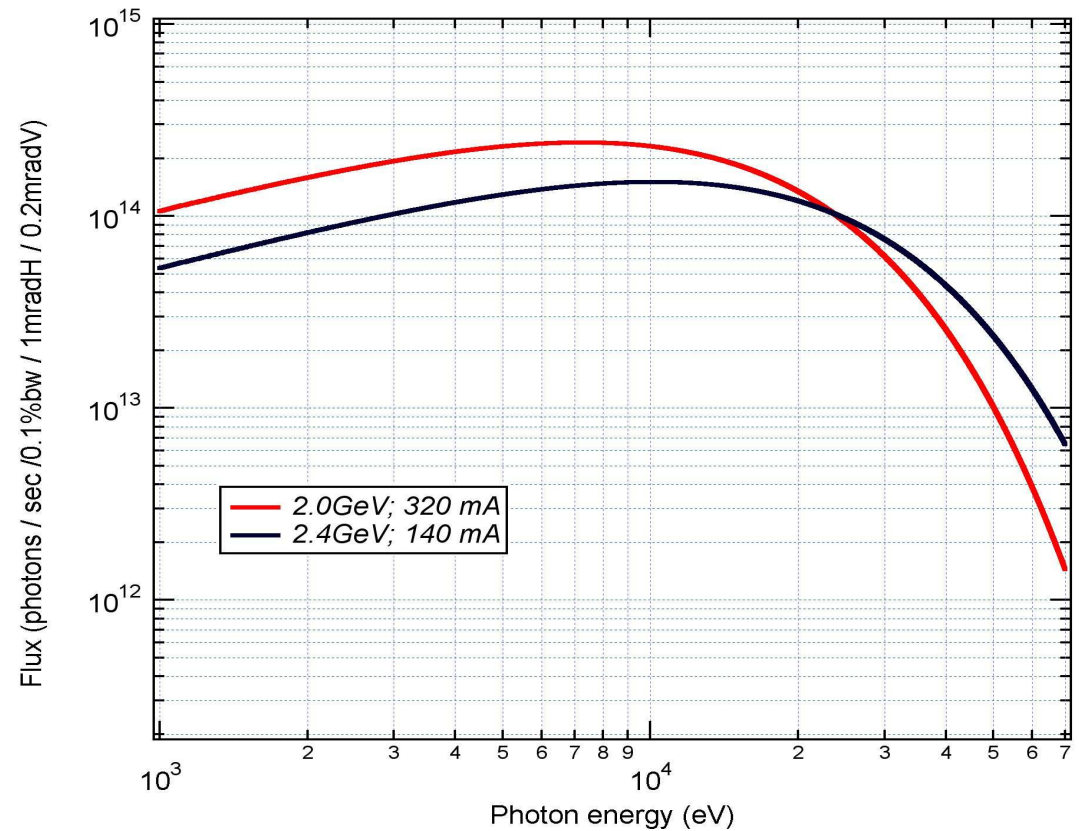

Sincrotrone Trieste S.c.p.A.
 Hard-X Ray Optics
 Laboratory

wavelength cutoff [Å]
 energy cutoff [keV]
 grazing angle[deg]

SCW 2GeV, 320 mA, 1 mradH x 0.2mradV

Incoming power [watt]
 absorbed power [watt]
 reflected power [watt]



mirror setup

	584	664	748
0.5 □ 24.8KeV 0.1505ū	180 404	185 479	190 557
0.55 □ 22.5KeV 0.167ū	213 371	219 446	225 523
0.6 □ 20.7KeV 0.185ū	250 334	257 407	265 483
0.73 □ 17KeV 0.225ū	329 255	339 325	349 398
	1.5 □ 8.3KeV 900μm	2 □ 6.2KeV 360μm	3 □ 4.1KeV 100μm




wavelength cutoff [Å]
 energy cutoff [keV]
 thickness [μm]

Graphite filter setup




Sincrotrone Trieste S.c.p.A.
 Hard-X Ray Optics
 Laboratory

ANSYS 5.7
NOV 20 2003
10:48:41


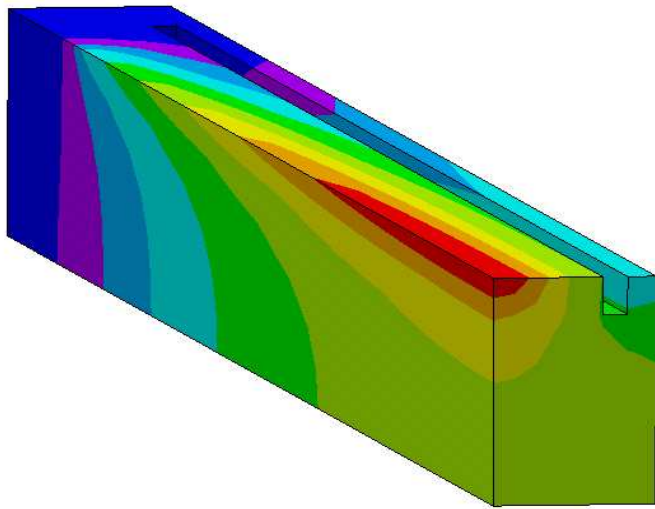


*case of C-filter 360 μ m / 0.185 $^\circ$ grazing angle
257 Watt absorbed*

Temperature distribution [K $^\circ$]

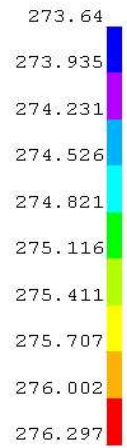


*Mirror deformation [m]
The coolant is water at 293K $^\circ$
and the flow is 300 l/h*

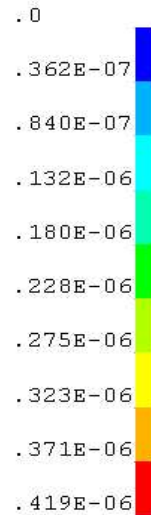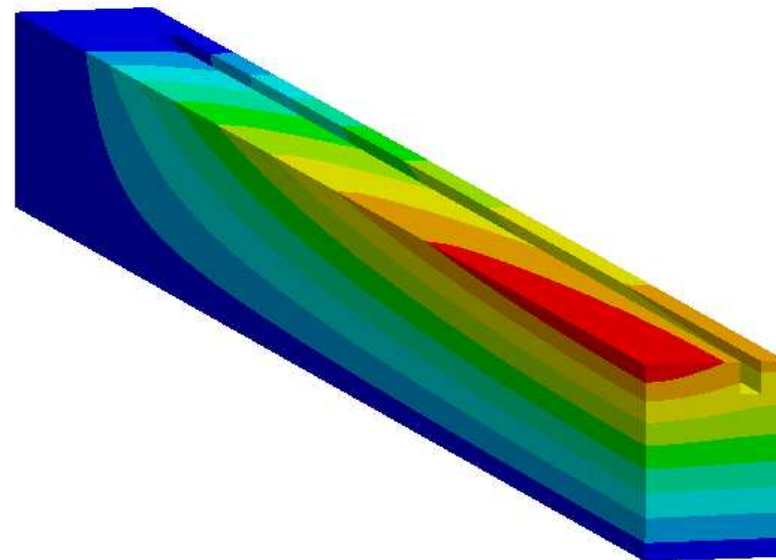

Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory


SCW 2GeV, 320 mA, 1 mradH x 0.2mradV

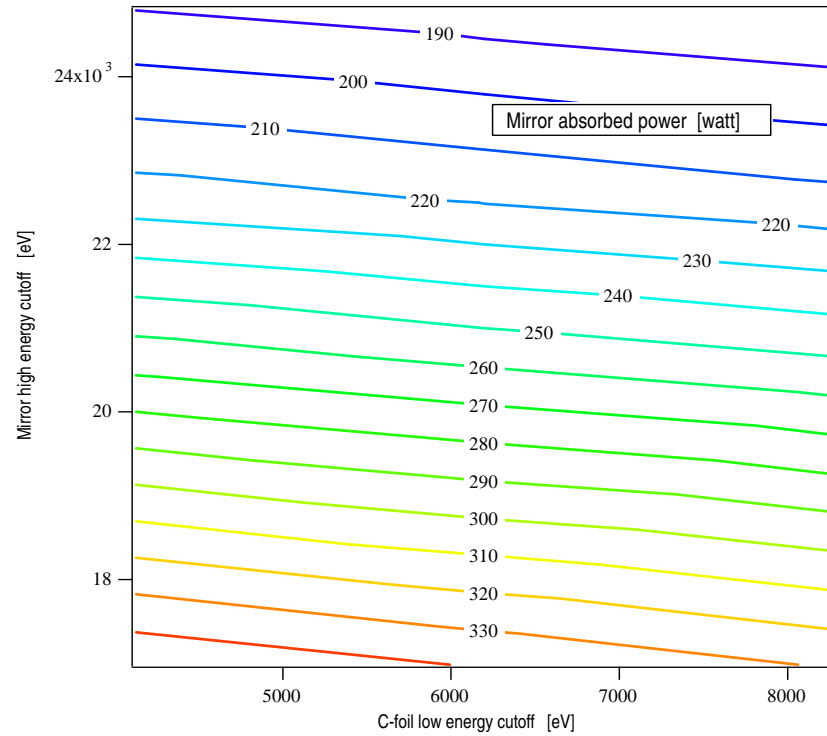


*because of the reflectivity the
the mirror absorbed power depends
on the low-pass filter threshold mostly*



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

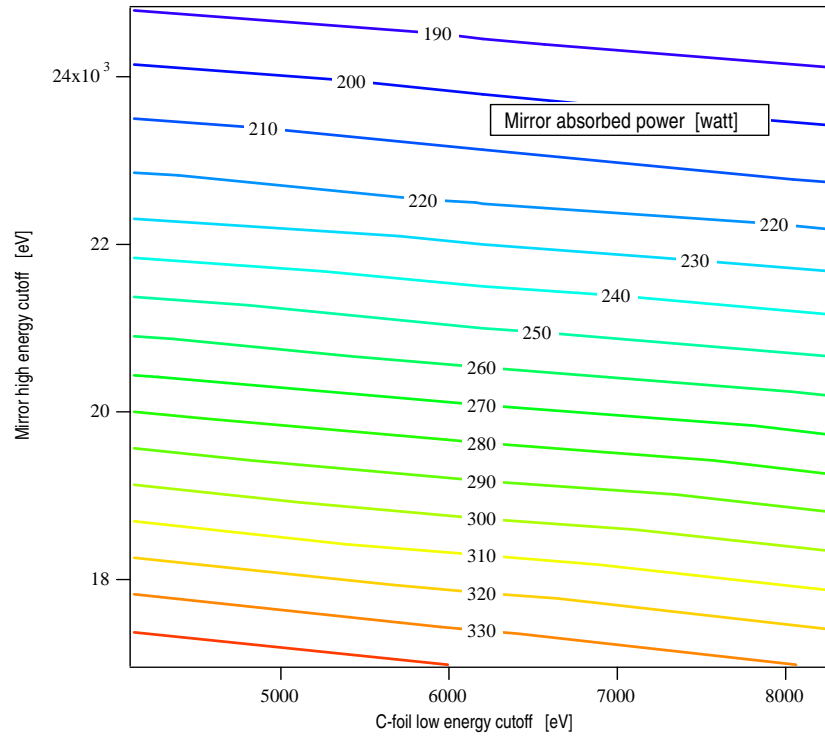

SCW 2GeV, 320 mA, 1 mradH x 0.2mradV

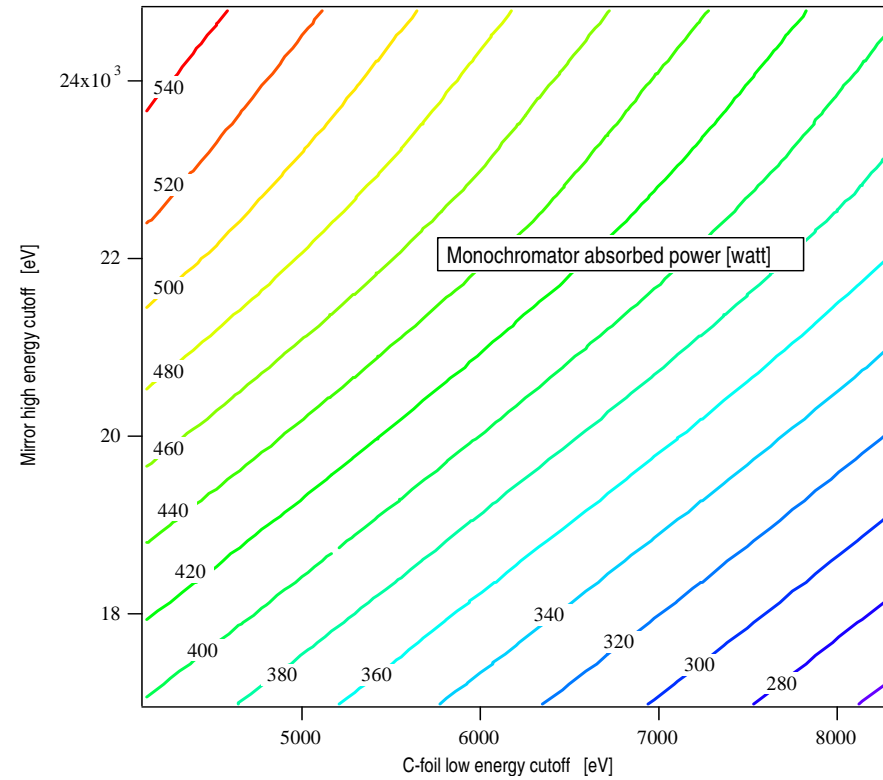


because of the reflectivity the mirror absorbed power depends on the low-pass filter threshold mostly

the high-pass filter threshold now has a not negligible contribution




Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

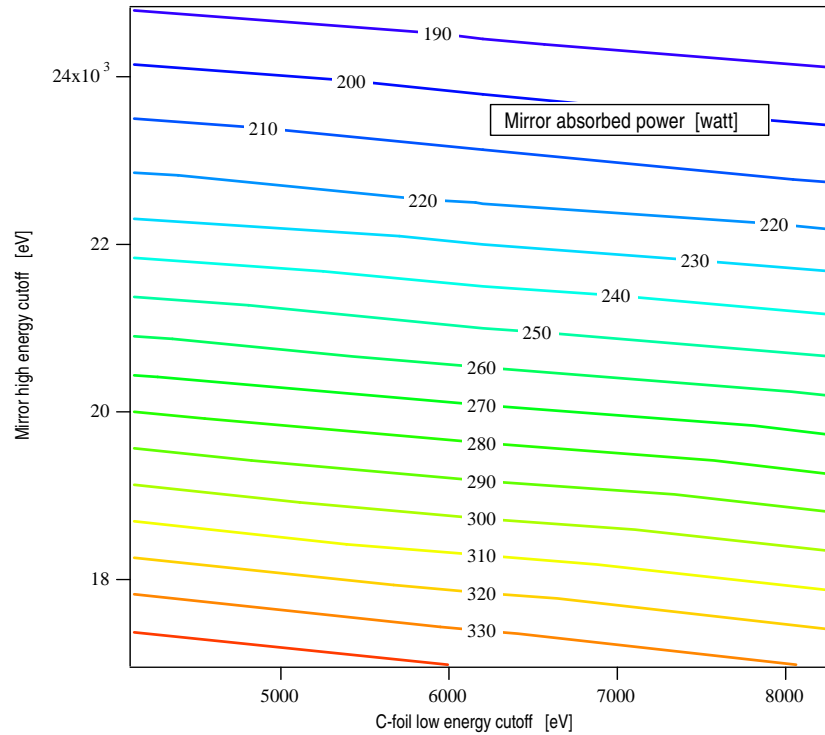

because of the reflectivity the mirror absorbed power depends on the low-pass filter threshold mostly



SCW 2GeV, 320 mA, 1 mradH x 0.2mradV

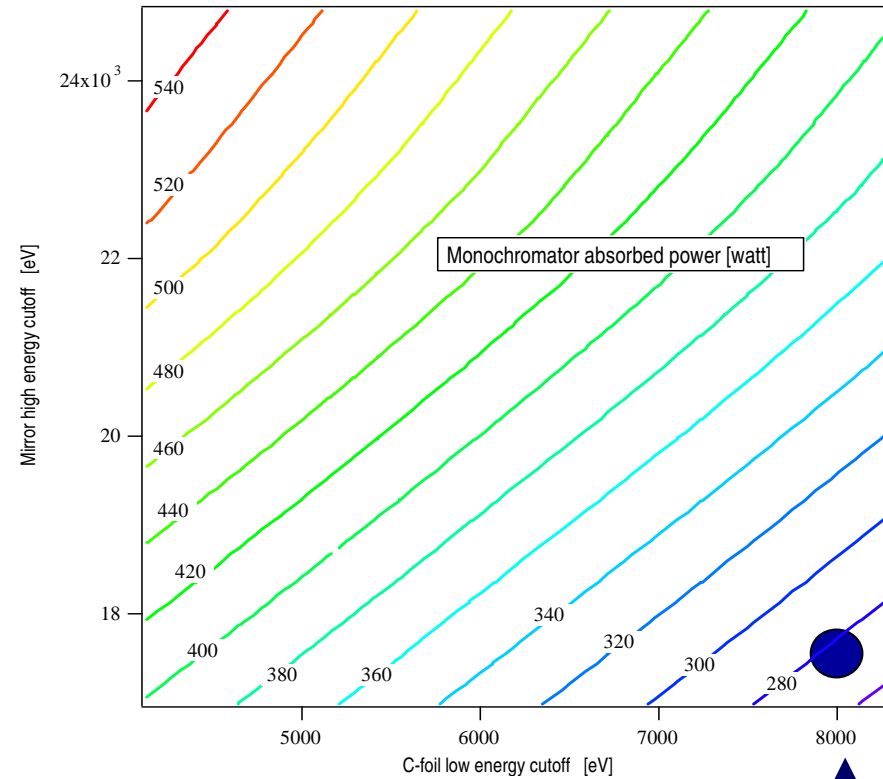


the high-pass filter threshold now has a not negligible contribution


first step: H₂O as coolant is still almost ok



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

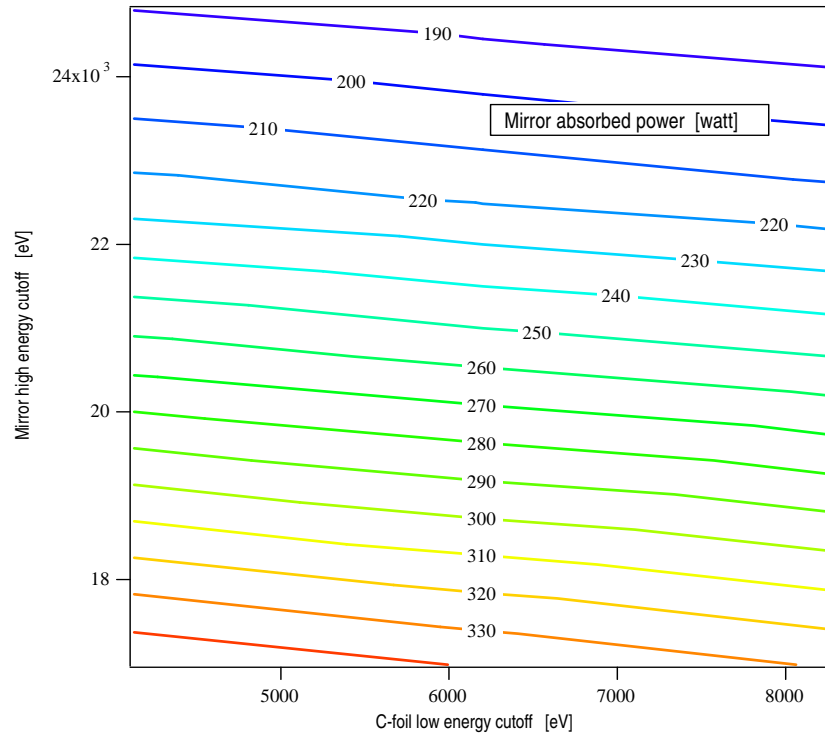

SCW 2GeV, 320 mA, 1 mradH x 0.2mradV

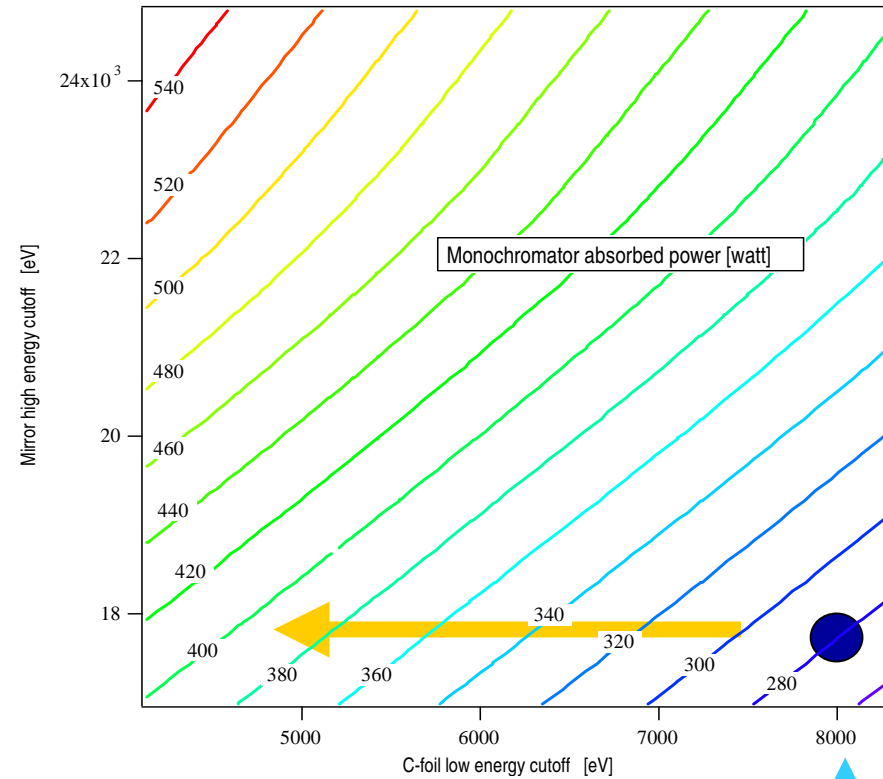

because of the reflectivity the mirror absorbed power depends on the low-pass filter threshold mostly


the high-pass filter threshold now has a not negligible contribution

first step: H₂O as coolant is still almost ok

in a second time: extend the range in the low energy direction




Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory

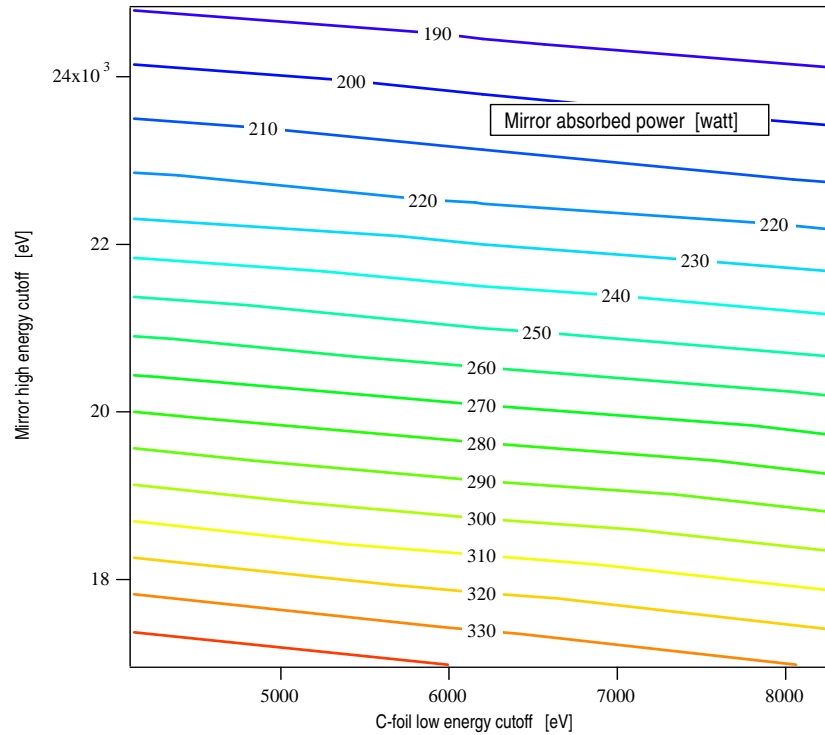

because of the reflectivity the mirror absorbed power depends on the low-pass filter threshold mostly



SCW 2GeV, 320 mA, 1 mradH x 0.2mradV

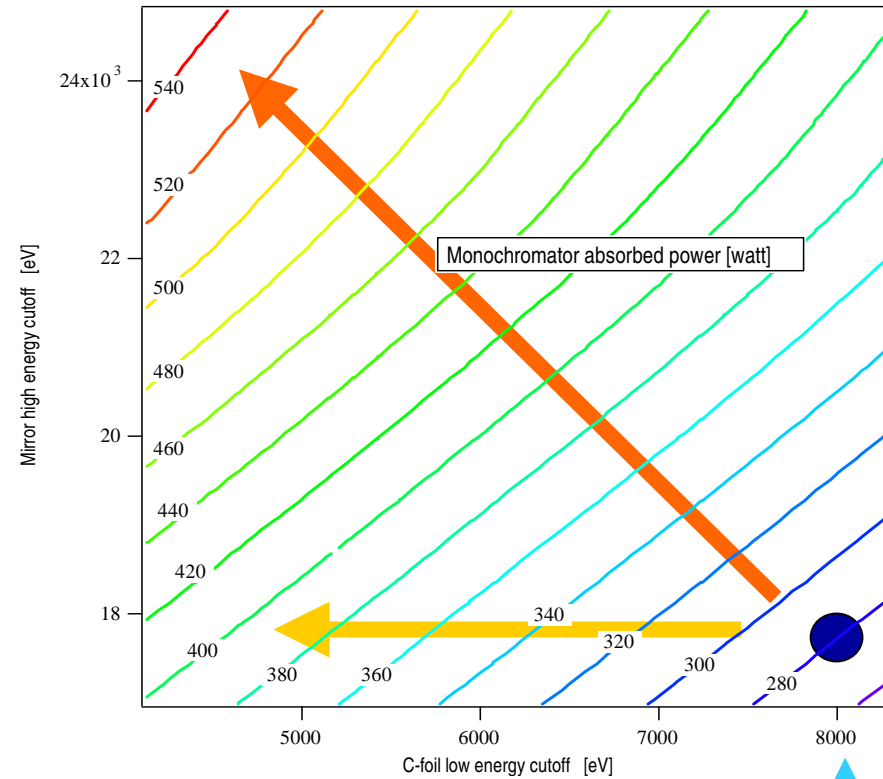

the high-pass filter threshold now has a not negligible contribution

first step: H₂O as coolant is still almost ok



in a second time: extend the range in the low energy direction or, if requested, open the energy window.....both with LN2 as coolant



Sincrotrone Trieste S.c.p.A.
Hard-X Ray Optics
Laboratory